

Anca-Fit/PINNACLE Total Conventional Hip Investigation

Note: This analysis compares the Anca-Fit/PINNACLE femoral stem/acetabular combination with all other total conventional hip prostheses.

This combination has been identified as having a significantly higher rate of revision. For a detailed explanation of the process used by the Registry that results in identification of prostheses that have a higher than anticipated rate of revision please refer to the Prostheses with Higher than Anticipated Rates of Revision chapter of the most recent AOANJRR Annual Report, <https://aoanjrr.sahmri.com/annual-reports-2026>.

Note: Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator. Procedures using prostheses with no recorded use in 2025 are excluded from the comparator.

TABLE 1

Revision Rate of Primary Total Conventional Hip Replacement

The revision rate of the Anca-Fit/PINNACLE total conventional hip combination is compared to all other total conventional hip prostheses.

Table 1: Revision Rates of Primary Total Conventional Hip Replacement

Component	N Revised	N Total	Obs. Years	Revisions/100 Obs. Yrs (95% CI)
Anca-Fit/PINNACLE	16	101	1152	1.39 (0.79, 2.25)
Other Total Conventional Hip	19701	568612	3623071	0.54 (0.54, 0.55)
TOTAL	19717	568713	3624223	0.54 (0.54, 0.55)

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 2

Yearly Cumulative Percent Revision of Primary Total Conventional Hip Replacement

The yearly cumulative percent revision of the Anca-Fit/PINNACLE total conventional hip combination is compared to all other total conventional hip prostheses.

Table 2: Yearly Cumulative Percent Revision (95% CI) of Primary Total Conventional Hip Replacement

CPR	1 Yr	2 Yrs	3 Yrs	4 Yrs	5 Yrs	6 Yrs
Anca-Fit/PINNACLE	6.0 (2.7, 12.8)	8.0 (4.1, 15.4)	8.0 (4.1, 15.4)	11.1 (6.3, 19.2)	11.1 (6.3, 19.2)	11.1 (6.3, 19.2)
Other Total Conventional Hip	1.7 (1.7, 1.8)	2.1 (2.1, 2.2)	2.5 (2.4, 2.5)	2.8 (2.7, 2.8)	3.0 (3.0, 3.1)	3.3 (3.2, 3.3)

CPR	7 Yrs	8 Yrs	9 Yrs	10 Yrs	11 Yrs	12 Yrs
Anca-Fit/PINNACLE	12.3 (7.2, 20.7)	13.5 (8.1, 22.2)	14.9 (9.1, 23.9)	16.4 (10.1, 25.8)	16.4 (10.1, 25.8)	16.4 (10.1, 25.8)
Other Total Conventional Hip	3.5 (3.5, 3.6)	3.8 (3.7, 3.9)	4.1 (4.0, 4.2)	4.4 (4.3, 4.5)	4.7 (4.6, 4.8)	5.1 (5.0, 5.2)

CPR	13 Yrs	14 Yrs	15 Yrs	16 Yrs	17 Yrs	18 Yrs
Anca-Fit/PINNACLE	16.4 (10.1, 25.8)	18.2 (11.4, 28.4)	18.2 (11.4, 28.4)			
Other Total Conventional Hip	5.4 (5.3, 5.5)	5.8 (5.7, 5.9)	6.2 (6.1, 6.3)	6.6 (6.5, 6.7)	7.0 (6.8, 7.1)	7.3 (7.1, 7.5)

CPR	19 Yrs	20 Yrs	21 Yrs	22 Yrs	23 Yrs	24 Yrs	25 Yrs
Anca-Fit/PINNACLE							
Other Total Conventional Hip	7.7 (7.5, 7.9)	8.0 (7.8, 8.2)	8.6 (8.3, 8.9)	9.0 (8.7, 9.4)	9.6 (9.1, 10.0)	9.6 (9.2, 10.1)	10.3 (9.4, 11.2)

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.


FIGURE 1

Yearly Cumulative Percent Revision of Primary Total Conventional Hip Replacement

The yearly cumulative percent revision of the Anca-Fit/PINNACLE total conventional hip combination is compared to all other total conventional hip prostheses. In addition, hazard ratios are reported.

Hazard ratios are reported for specific time periods during which the hazard ratio is constant. This is done to enable more specific and valid comparisons of the risk of revision over time. The pattern of variation in risk has important implications with respect to the underlying reasons for any difference.

Figure 1: Cumulative Percent Revision of Primary Total Conventional Hip Replacement

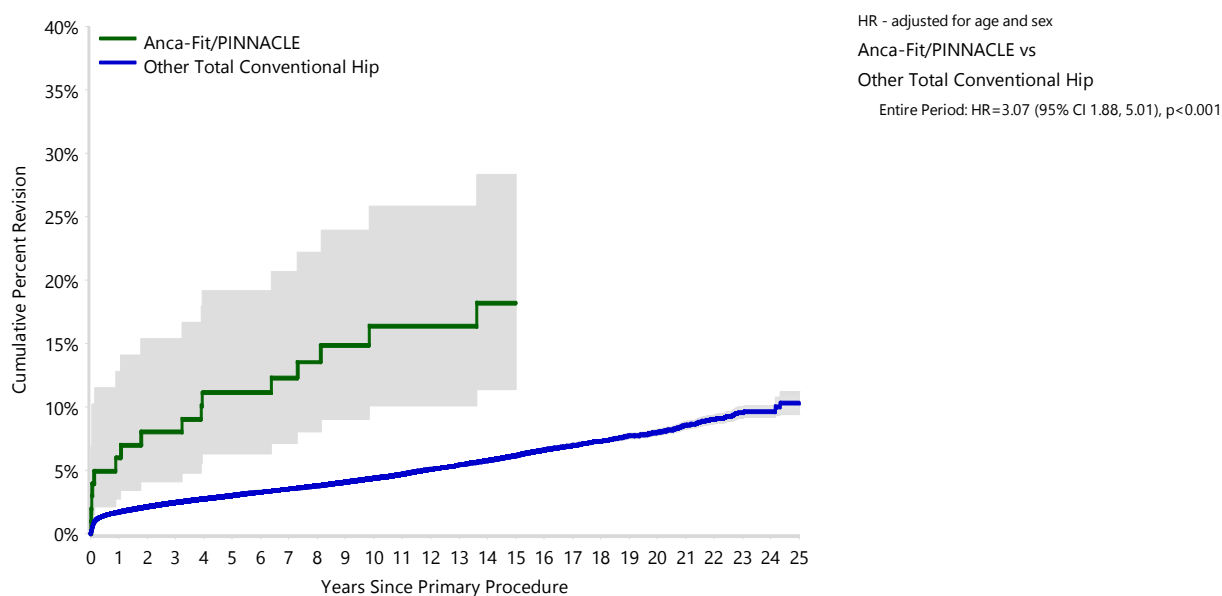

Number at Risk	0 Yr	1 Yr	2 Yrs	3 Yrs	4 Yrs	5 Yrs	6 Yrs	7 Yrs
Anca-Fit/PINNACLE	101	92	90	90	86	80	77	72
Other Total Conventional Hip	568612	501603	445433	391564	344389	299354	259516	220948

Number at Risk	8 Yrs	9 Yrs	10 Yrs	11 Yrs	12 Yrs	13 Yrs	14 Yrs	15 Yrs	16 Yrs
Anca-Fit/PINNACLE	65	60	57	51	50	50	44	40	36
Other Total Conventional Hip	185798	153739	124848	99575	79802	63343	49730	38033	28355

Number at Risk	17 Yrs	18 Yrs	19 Yrs	20 Yrs	21 Yrs	22 Yrs	23 Yrs	24 Yrs	25 Yrs
Anca-Fit/PINNACLE	34	27	6	0	0	0	0	0	0
Other Total Conventional Hip	20437	14638	10400	7396	4859	2863	1455	560	95

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 3**Primary Diagnosis for Revised Primary Total Conventional Hip Replacement**

This table identifies the diagnosis of the primary procedure which was subsequently revised. This information is provided as there is a variation on outcome depending on the primary diagnosis. It is therefore important when considering the reasons for a higher than anticipated rate of revision that there is identification of the primary diagnosis. This information should be compared to the primary diagnosis for the revisions of all other total conventional hip prostheses.

Table 3: Primary Diagnosis for Revised Primary Total Conventional Hip Replacement

Primary Diagnosis	Anca-Fit/PINNACLE		Other Total Conventional Hip	
	Number	Percent	Number	Percent
Osteoarthritis	11	68.8	16296	82.7
Fractured Neck Of Femur	1	6.3	1521	7.7
Osteonecrosis	4	25.0	844	4.3
Developmental Dysplasia			321	1.6
Rheumatoid Arthritis			201	1.0
Failed Internal Fixation			165	0.8
Tumour			157	0.8
Other Inflammatory Arthritis			102	0.5
Fracture/Dislocation			60	0.3
Other			21	0.1
Arthrodesis Takedown			13	0.1
TOTAL	16	100.0	19701	100.0

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 4

Reasons for Revision

This is reported in two ways: a percentage of primary procedures revised and as a percentage of all revision procedures.

% Primaries Revised: This shows the proportional contribution of each revision diagnosis as a percentage of the total number of primary procedures. This percentage can be used to approximate the risk of being revised for that diagnosis. Differing percentages between groups, with the same distribution of follow up time, may identify problems of concern.

% Revisions: The number of revisions for each diagnosis is expressed as a percentage of the total number of revisions. This shows the distribution of reasons for revision within a group but cannot be used as a comparison between groups.

Table 4: Primary Total Conventional Hip Replacement - Reason for Revision (Follow-up Limited to 19.6 Years)

Revision Diagnosis	Anca-Fit/PINNACLE			Other Total Conventional Hip		
	Number	% Primaries Revised	% Revisions	Number	% Primaries Revised	% Revisions
Infection	1	1.0	6.3	5005	0.9	25.5
Prosthesis Dislocation/Instability	2	2.0	12.5	4374	0.8	22.3
Fracture	9	8.9	56.3	4334	0.8	22.1
Loosening	4	4.0	25.0	3606	0.6	18.4
Pain				339	0.1	1.7
Leg Length Discrepancy				307	0.1	1.6
Malposition				281	0.0	1.4
Implant Breakage Stem				208	0.0	1.1
Lysis				194	0.0	1.0
Implant Breakage Acetabular Insert				126	0.0	0.6
Wear Acetabular Insert				103	0.0	0.5
Incorrect Sizing				96	0.0	0.5
Metal Related Pathology				89	0.0	0.5
Implant Breakage Acetabular				73	0.0	0.4
Wear Head				48	0.0	0.2
Tumour				46	0.0	0.2
Heterotopic Bone				26	0.0	0.1
Implant Breakage Head				26	0.0	0.1
Wear Acetabulum				6	0.0	0.0
Osteonecrosis				4	0.0	0.0
Synovitis				2	0.0	0.0
Progression Of Disease				1	0.0	0.0
Other				321	0.1	1.6
N Revision	16	15.8	100.0	19615	3.4	100.0
N Primary	101			568612		

Note: This table is restricted to revisions within 19.6 years for all groups to allow a time-matched comparison of revisions.

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

FIGURE 2**Cumulative Incidence Revision Diagnosis of Primary Total Conventional Hip Replacement**

This figure details the cumulative incidence of the most common reasons for revision. The five most common reasons for revision are included as long as each of these reasons account for more than 10 procedures or at least 5% of all revisions for the Anca-Fit/PINNACLE total conventional hip combination. A comparative graph is provided of the cumulative incidence for the same reasons for revisions for all other total conventional hip prostheses.

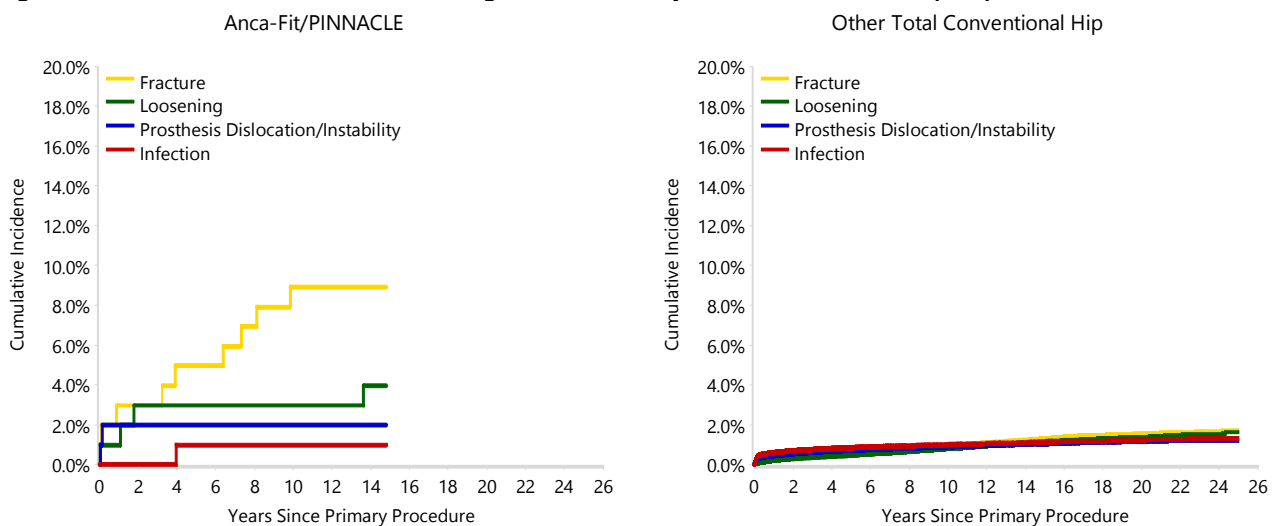
Figure 2: Cumulative Incidence Revision Diagnosis for Primary Total Conventional Hip Replacement

TABLE 5

Type of Revision Performed for Primary Total Conventional Hip Replacement

This analysis identifies the components used in the revision of the Anca-Fit/PINNACLE total conventional hip combination and compares it to the components used in the revision of all other total conventional hip prostheses.

The reason this analysis is undertaken is to identify whether there is one or more components which are being replaced that differ from the components replaced for revisions of all other total conventional hip prostheses i.e. is there a difference in the type of revision undertaken for the Anca-Fit/PINNACLE total conventional hip combination compared to all other total conventional hip prostheses.

Table 5: Primary Total Conventional Hip Replacement - Type of Revision (Follow-up Limited to 19.6 Years)

Type of Revision	Anca-Fit/PINNACLE		Other Total Conventional Hip	
	Number	Percent	Number	Percent
Femoral Component	10	62.5	6627	33.8
Acetabular Component			3293	16.8
THR (Femoral/Acetabular)	1	6.3	2280	11.6
Cement Spacer			556	2.8
Removal of Prostheses			92	0.5
Reinsertion of Components			27	0.1
Total Femoral			12	0.1
Bipolar Head and Femoral			10	0.1
N Major	11	68.8	12897	65.8
Head/Insert	1	6.3	5351	27.3
Head Only			911	4.6
Minor Components	1	6.3	285	1.5
Insert Only			168	0.9
Head/Neck	3	18.8		
Bipolar Only			1	0.0
Cement Only			1	0.0
Head/Neck/Insert			1	0.0
N Minor	5	31.3	6718	34.2
TOTAL	16	100.0	19615	100.0

Note: This table is restricted to revisions within 19.6 years for all groups to allow a time-matched comparison of revisions.

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 6**Revision Rates of Anca-Fit/PINNACLE Primary Total Conventional Hip Replacement by Fixation**

This analysis is provided as some prostheses have more than one fixation option. Additionally there are prostheses where an alternative to the recommended approach to fixation was used e.g. a cementless prosthesis that has been cemented or vice-versa.

Table 6: Revised Number of Anca-Fit/PINNACLE Primary Total Conventional Hip Replacement by Fixation

Fixation	N Revised	N Total
Cementless	16	101
TOTAL	16	101

TABLE 7**Revision Rates of Anca-Fit/PINNACLE Primary Total Conventional Hip Replacement by Bearing Surface**

This analysis is provided as some prostheses are combined with a variety of bearing surfaces. All bearing surfaces used with this combination are listed.

Table 7: Revised Number of Anca-Fit/PINNACLE Primary Total Conventional Hip Replacement by Bearing Surface

Bearing Surface	N Revised	N Total
Ceramic/Ceramic	15	98
Ceramic/XLPE	0	1
Metal/XLPE	1	2
TOTAL	16	101

TABLE 8

Revision Rates of Primary Total Conventional Hip Replacement by State

This enables a state by state variation to be identified for the Anca-Fit/PINNACLE total conventional hip combination and provides the comparative data for each of the states for all other total conventional hip prostheses.

The purpose of this analysis is to determine if the higher than anticipated rate of revision has widespread distribution between states. If there is widespread distribution then the reason for the higher than anticipated rate of revision is unlikely to be surgeon specific. If the prosthesis has been used in only a small number of states it is not possible to distinguish if the higher than anticipated rate of revision is related to the prosthesis, surgeon, technique or patient.

Table 8: Revised Number of Primary Total Conventional Hip Replacement by State

Component	State	N Revised	N Total
Anca-Fit/PINNACLE	NSW	16	101
Other Total Conventional Hip	NSW	5507	167002
	VIC	4853	146558
	QLD	3968	104359
	WA	2510	64946
	SA	1887	52084
	TAS	439	18269
	ACT/NT	537	15394
TOTAL		19717	568713

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 9**Number of Revisions of Anca-Fit/PINNACLE Primary Total Conventional Hip Replacement by Year of Implant**

This analysis details the number of prostheses reported each year to the Registry for the Anca-Fit/PINNACLE total conventional hip combination. It also provides the subsequent number of revisions of the primaries reported in that year.

Primary procedures performed in later years have had less follow up time therefore the number revised is expected to be less than the number revised in earlier years. For example, a primary procedure performed in 2025 has a maximum of one year to be revised, whereas a primary procedure performed in 2023 has a maximum of three years to be revised.

Table 9: Number of Revisions of Anca-Fit/PINNACLE Primary Total Conventional Hip Replacement by Year of Implant

Year of Implant	Number Revised	Total Number
2006	6	30
2007	7	55
2008	3	16
TOTAL	16	101

TABLE 10

Revision Rates of Anca-Fit/PINNACLE Primary Total Conventional Hip Replacement by Catalogue Number Range

Many prostheses have a number of catalogue ranges. The catalogue range is specific to particular design features; more than one catalogue range usually indicates a minor difference in design in a particular Anca-Fit/PINNACLE prosthesis.

This analysis has been undertaken to determine if the revision rate varies according to the catalogue number range.

Model	Catalogue Range	Catalogue Description	Cement	Material	Stem Type	Coating	Fixation
Femoral Stem							
Anca-Fit	PPR67600-PPR67630	1/3 HA GRT BLAST ANATOMICAL STEM	NO	METAL	REQUIRES FEMNECK	HA COATED	
Acetabular							
PINNACLE	121712048-121712066	SECTOR ACETABULAR CUP W/DUOFIX HA	NO	METAL		HA COATED	
PINNACLE	121720048-121720072	MULTIHOLE ACETABULAR SHELL W/POROCOAT	NO	METAL			POROUS

Table 10: Revised Number of Anca-Fit/PINNACLE Primary Total Conventional Hip Replacement by Catalogue Number Range

Femoral Stem Range	Acetabular Range	N Revised	N Total
PPR67600-PPR67630	121712048-121712066	16	100
	121720048-121720072	0	1
TOTAL		16	101