

Novation Total Conventional Hip Investigation

Note: This analysis compares the Novation femoral stem prosthesis with all other total conventional hip prostheses.

This prosthesis has been identified as having a significantly higher rate of revision. For a detailed explanation of the process used by the Registry that results in identification of prostheses that have a higher than anticipated rate of revision please refer to the Prostheses with Higher than Anticipated Rates of Revision chapter of the most recent AOANJRR Annual Report, <https://aoanjrr.sahmri.com/annual-reports-2026>.

Note: Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator. Procedures using prostheses with no recorded use in 2025 are excluded from the comparator.

TABLE 1

Revision Rate of Primary Total Conventional Hip Replacement

The revision rate of the Novation total conventional hip prosthesis is compared to all other total conventional hip prostheses.

Table 1: Revision Rates of Primary Total Conventional Hip Replacement

Component	N Revised	N Total	Obs. Years	Revisions/100 Obs. Yrs (95% CI)
Novation	100	1423	12715	0.79 (0.64, 0.96)
Other Total Conventional Hip	19701	568612	3623071	0.54 (0.54, 0.55)
TOTAL	19801	570035	3635786	0.54 (0.54, 0.55)

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 2

Yearly Cumulative Percent Revision of Primary Total Conventional Hip Replacement

The yearly cumulative percent revision of the Novation total conventional hip prosthesis is compared to all other total conventional hip prostheses.

Table 2: Yearly Cumulative Percent Revision (95% CI) of Primary Total Conventional Hip Replacement

CPR	1 Yr	2 Yrs	3 Yrs	4 Yrs	5 Yrs	6 Yrs
Novation	3.7 (2.9, 4.9)	4.2 (3.2, 5.3)	4.4 (3.5, 5.7)	4.9 (3.9, 6.1)	5.2 (4.1, 6.5)	5.4 (4.4, 6.7)
Other Total Conventional Hip	1.7 (1.7, 1.8)	2.1 (2.1, 2.2)	2.5 (2.4, 2.5)	2.8 (2.7, 2.8)	3.0 (3.0, 3.1)	3.3 (3.2, 3.3)

CPR	7 Yrs	8 Yrs	9 Yrs	10 Yrs	11 Yrs	12 Yrs
Novation	5.8 (4.7, 7.2)	6.3 (5.2, 7.8)	7.1 (5.8, 8.7)	7.5 (6.1, 9.1)	7.8 (6.4, 9.5)	8.3 (6.8, 10.1)
Other Total Conventional Hip	3.5 (3.5, 3.6)	3.8 (3.7, 3.9)	4.1 (4.0, 4.2)	4.4 (4.3, 4.5)	4.7 (4.6, 4.8)	5.1 (5.0, 5.2)

CPR	13 Yrs	14 Yrs	15 Yrs	16 Yrs	17 Yrs	18 Yrs
Novation	8.3 (6.8, 10.1)	8.3 (6.8, 10.1)				
Other Total Conventional Hip	5.4 (5.3, 5.5)	5.8 (5.7, 5.9)	6.2 (6.1, 6.3)	6.6 (6.5, 6.7)	7.0 (6.8, 7.1)	7.3 (7.1, 7.5)

CPR	19 Yrs	20 Yrs	21 Yrs	22 Yrs	23 Yrs	24 Yrs	25 Yrs
Novation							
Other Total Conventional Hip	7.7 (7.5, 7.9)	8.0 (7.8, 8.2)	8.6 (8.3, 8.9)	9.0 (8.7, 9.4)	9.6 (9.1, 10.0)	9.6 (9.2, 10.1)	10.3 (9.4, 11.2)

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.


FIGURE 1

Yearly Cumulative Percent Revision of Primary Total Conventional Hip Replacement

The yearly cumulative percent revision of the Novation total conventional hip prosthesis is compared to all other total conventional hip prostheses. In addition, hazard ratios are reported.

Hazard ratios are reported for specific time periods during which the hazard ratio is constant. This is done to enable more specific and valid comparisons of the risk of revision over time. The pattern of variation in risk has important implications with respect to the underlying reasons for any difference.

Figure 1: Cumulative Percent Revision of Primary Total Conventional Hip Replacement

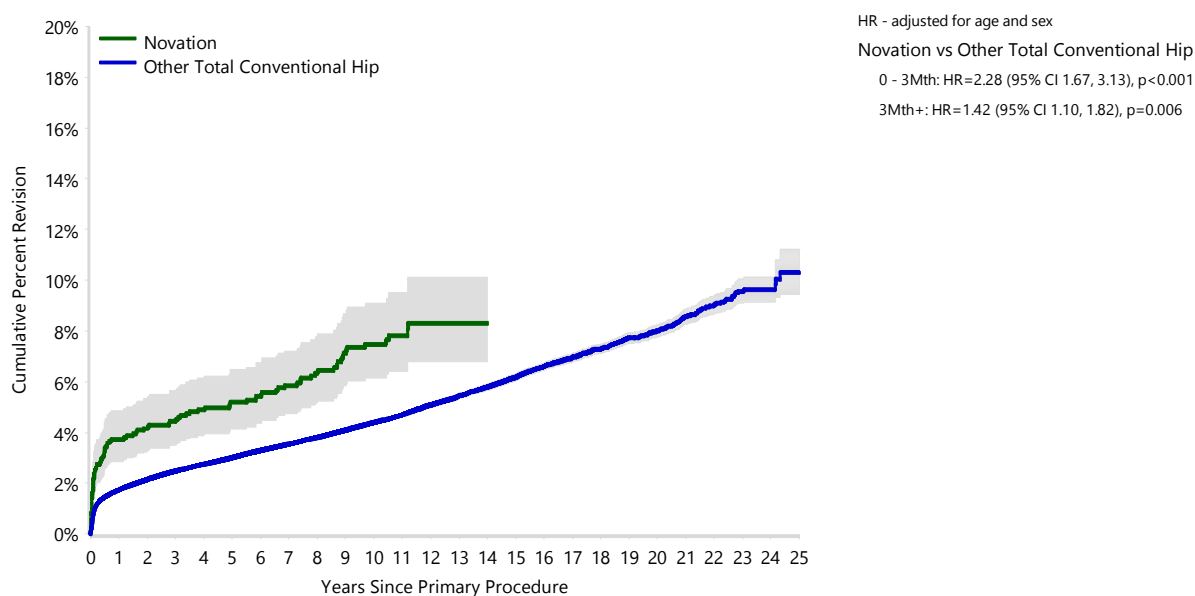

Number at Risk	0 Yr	1 Yr	2 Yrs	3 Yrs	4 Yrs	5 Yrs	6 Yrs	7 Yrs
Novation	1423	1360	1343	1319	1285	1240	1176	1025
Other Total Conventional Hip	568612	501603	445433	391564	344389	299354	259516	220948

Number at Risk	8 Yrs	9 Yrs	10 Yrs	11 Yrs	12 Yrs	13 Yrs	14 Yrs	15 Yrs	16 Yrs
Novation	908	803	658	422	243	144	57	23	3
Other Total Conventional Hip	185798	153739	124848	99575	79802	63343	49730	38033	28355

Number at Risk	17 Yrs	18 Yrs	19 Yrs	20 Yrs	21 Yrs	22 Yrs	23 Yrs	24 Yrs	25 Yrs
Novation	0	0	0	0	0	0	0	0	0
Other Total Conventional Hip	20437	14638	10400	7396	4859	2863	1455	560	95

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 3**Primary Diagnosis for Revised Primary Total Conventional Hip Replacement**

This table identifies the diagnosis of the primary procedure which was subsequently revised. This information is provided as there is a variation on outcome depending on the primary diagnosis. It is therefore important when considering the reasons for a higher than anticipated rate of revision that there is identification of the primary diagnosis. This information should be compared to the primary diagnosis for the revisions of all other total conventional hip prostheses.

Table 3: Primary Diagnosis for Revised Primary Total Conventional Hip Replacement

Primary Diagnosis	Novation		Other Total Conventional Hip	
	Number	Percent	Number	Percent
Osteoarthritis	92	92.0	16296	82.7
Fractured Neck Of Femur	3	3.0	1521	7.7
Osteonecrosis	4	4.0	844	4.3
Developmental Dysplasia			321	1.6
Rheumatoid Arthritis			201	1.0
Failed Internal Fixation			165	0.8
Tumour	1	1.0	157	0.8
Other Inflammatory Arthritis			102	0.5
Fracture/Dislocation			60	0.3
Other			21	0.1
Arthrodesis Takedown			13	0.1
TOTAL	100	100.0	19701	100.0

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 4

Reasons for Revision

This is reported in two ways: a percentage of primary procedures revised and as a percentage of all revision procedures.

% Primaries Revised: This shows the proportional contribution of each revision diagnosis as a percentage of the total number of primary procedures. This percentage can be used to approximate the risk of being revised for that diagnosis. Differing percentages between groups, with the same distribution of follow up time, may identify problems of concern.

% Revisions: The number of revisions for each diagnosis is expressed as a percentage of the total number of revisions. This shows the distribution of reasons for revision within a group but cannot be used as a comparison between groups.

Table 4: Primary Total Conventional Hip Replacement - Reason for Revision (Follow-up Limited to 16.2 Years)

Revision Diagnosis	Number	Novation		Other Total Conventional Hip		
		% Primaries Revised	% Revisions	Number	% Primaries Revised	% Revisions
Infection	13	0.9	13.0	4983	0.9	25.7
Prosthesis Dislocation/Instability	19	1.3	19.0	4347	0.8	22.4
Fracture	29	2.0	29.0	4286	0.8	22.1
Loosening	27	1.9	27.0	3552	0.6	18.3
Pain				337	0.1	1.7
Leg Length Discrepancy	2	0.1	2.0	307	0.1	1.6
Malposition	4	0.3	4.0	280	0.0	1.4
Implant Breakage Stem	1	0.1	1.0	202	0.0	1.0
Lysis	1	0.1	1.0	180	0.0	0.9
Implant Breakage Acetabular Insert	1	0.1	1.0	125	0.0	0.6
Incorrect Sizing	1	0.1	1.0	96	0.0	0.5
Wear Acetabular Insert	1	0.1	1.0	94	0.0	0.5
Metal Related Pathology				80	0.0	0.4
Implant Breakage Acetabular				71	0.0	0.4
Wear Head				47	0.0	0.2
Tumour	1	0.1	1.0	46	0.0	0.2
Heterotopic Bone				26	0.0	0.1
Implant Breakage Head				23	0.0	0.1
Wear Acetabulum				6	0.0	0.0
Osteonecrosis				4	0.0	0.0
Synovitis				2	0.0	0.0
Progression Of Disease				1	0.0	0.0
Other				320	0.1	1.6
N Revision	100	7.0	100.0	19415	3.4	100.0
N Primary	1423			568612		

Note: This table is restricted to revisions within 16.2 years for all groups to allow a time-matched comparison of revisions.

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

FIGURE 2**Cumulative Incidence Revision Diagnosis of Primary Total Conventional Hip Replacement**

This figure details the cumulative incidence of the most common reasons for revision. The five most common reasons for revision are included as long as each of these reasons account for more than 10 procedures or at least 5% of all revisions for the Novation total conventional hip prosthesis. A comparative graph is provided of the cumulative incidence for the same reasons for revisions for all other total conventional hip prostheses.

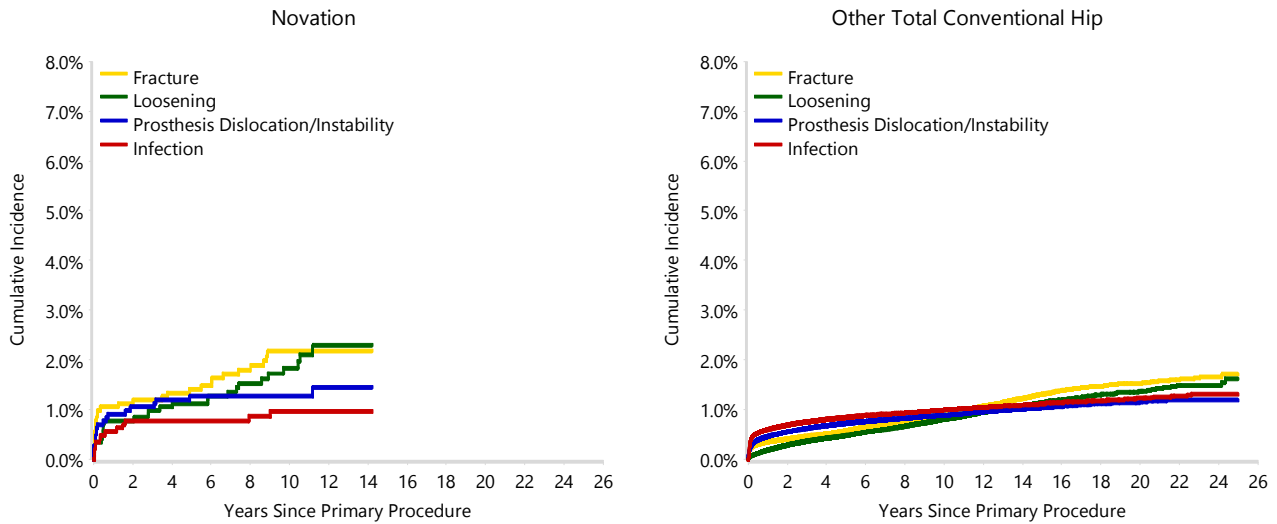
Figure 2: Cumulative Incidence Revision Diagnosis for Primary Total Conventional Hip Replacement

TABLE 5

Type of Revision Performed for Primary Total Conventional Hip Replacement

This analysis identifies the components used in the revision of the Novation total conventional hip prosthesis and compares it to the components used in the revision of all other total conventional hip prostheses.

The reason this analysis is undertaken is to identify whether there is one or more components which are being replaced that differ from the components replaced for revisions of all other total conventional hip prostheses i.e. is there a difference in the type of revision undertaken for the Novation total conventional hip prosthesis compared to all other total conventional hip prostheses.

Table 5: Primary Total Conventional Hip Replacement - Type of Revision (Follow-up Limited to 16.2 Years)

Type of Revision	Novation		Other Total Conventional Hip	
	Number	Percent	Number	Percent
Femoral Component	38	38.0	6565	33.8
Acetabular Component	24	24.0	3238	16.7
THR (Femoral/Acetabular)	11	11.0	2240	11.5
Cement Spacer	4	4.0	555	2.9
Removal of Prostheses			91	0.5
Reinsertion of Components			27	0.1
Total Femoral			12	0.1
Bipolar Head and Femoral			10	0.1
N Major	77	77.0	12738	65.6
Head/Insert	14	14.0	5319	27.4
Head Only	7	7.0	907	4.7
Minor Components	2	2.0	281	1.4
Insert Only			167	0.9
Bipolar Only			1	0.0
Cement Only			1	0.0
Head/Neck/Insert			1	0.0
N Minor	23	23.0	6677	34.4
TOTAL	100	100.0	19415	100.0

Note: This table is restricted to revisions within 16.2 years for all groups to allow a time-matched comparison of revisions.

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 6**Revision Rates of Novation Primary Total Conventional Hip Replacement by Fixation**

This analysis is provided as some prostheses have more than one fixation option. Additionally there are prostheses where an alternative to the recommended approach to fixation was used e.g. a cementless prosthesis that has been cemented or vice-versa.

Table 6: Revised Number of Novation Primary Total Conventional Hip Replacement by Fixation

Fixation	N Revised	N Total
Cementless	98	1417
Hybrid (Femur Cemented)	1	3
Reverse Hybrid (Femur Cementless)	1	3
TOTAL	100	1423

TABLE 7**Revision Rates of Novation Primary Total Conventional Hip Replacement by Bearing Surface**

This analysis is provided as some prostheses are combined with a variety of bearing surfaces. All bearing surfaces used with this prosthesis are listed.

Table 7: Revised Number of Novation Primary Total Conventional Hip Replacement by Bearing Surface

Bearing Surface	N Revised	N Total
Ceramic/Ceramic	42	637
Ceramic/Non XLPE	16	169
Ceramic/XLPE	30	477
Ceramic/XLPE + Antioxidant	0	2
Metal/Non XLPE	3	26
Metal/XLPE	9	110
Metal/XLPE + Antioxidant	0	1
Unknown	0	1
TOTAL	100	1423

TABLE 8**Revision Rates of Novation Primary Total Conventional Hip Replacement by Approach**

This analysis is provided as some prostheses are used with a variety of surgical approaches. All surgical approaches used with this prosthesis are listed.

Table 8: Revised Number of Novation Primary Total Conventional Hip Replacement by Approach

Approach	N Revised	N Total
Anterior	19	270
Lateral	9	159
Posterior	22	341
TOTAL	50	770

Note: Excludes 653 procedures with no approach recorded

TABLE 9

Revision Rates of Primary Total Conventional Hip Replacement by State

This enables a state by state variation to be identified for the Novation total conventional hip prosthesis and provides the comparative data for each of the states for all other total conventional hip prostheses.

The purpose of this analysis is to determine if the higher than anticipated rate of revision has widespread distribution between states. If there is widespread distribution then the reason for the higher than anticipated rate of revision is unlikely to be surgeon specific. If the prosthesis has been used in only a small number of states it is not possible to distinguish if the higher than anticipated rate of revision is related to the prosthesis, surgeon, technique or patient.

Table 9: Revised Number of Primary Total Conventional Hip Replacement by State

Component	State	N Revised	N Total
Novation	NSW	38	530
	VIC	14	84
	QLD	13	164
	WA	26	396
	TAS	9	249
Other Total Conventional Hip	NSW	5507	167002
	VIC	4853	146558
	QLD	3968	104359
	WA	2510	64946
	SA	1887	52084
	TAS	439	18269
	ACT/NT	537	15394
TOTAL		19801	570035

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 10**Number of Revisions of Novation Primary Total Conventional Hip Replacement by Year of Implant**

This analysis details the number of prostheses reported each year to the Registry for the Novation total conventional hip prosthesis. It also provides the subsequent number of revisions of the primaries reported in that year.

Primary procedures performed in later years have had less follow up time therefore the number revised is expected to be less than the number revised in earlier years. For example, a primary procedure performed in 2025 has a maximum of one year to be revised, whereas a primary procedure performed in 2023 has a maximum of three years to be revised.

Table 10: Number of Revisions of Novation Primary Total Conventional Hip Replacement by Year of Implant

Year of Implant	Number Revised	Total Number
2009	0	4
2010	1	32
2011	7	53
2012	13	130
2013	10	137
2014	13	226
2015	15	266
2016	8	148
2017	6	90
2018	9	101
2019	14	145
2020	3	51
2021	1	19
2022	0	21
TOTAL	100	1423

TABLE 11

Revision Rates of Novation Primary Total Conventional Hip Replacement by Catalogue Number Range

Many prostheses have a number of catalogue ranges. The catalogue range is specific to particular design features; more than one catalogue range usually indicates a minor difference in design in a particular Novation prosthesis.

This analysis has been undertaken to determine if the revision rate varies according to the catalogue number range.

Model	Catalogue Range	Catalogue Description	Cement	Material	Coating
Femoral Stem					
Novation	1600009-1600018	PLASMA TAPERED PRESS-FIT STANDARD OFFSET FEMORAL STEM	NO	METAL	HA COATED
Novation	1600109-1600118	PLASMA TAPERED PRESS-FIT EXTENDED OFFSET FEMORAL STEM	NO	METAL	HA COATED
Novation	1640108-1640118	HA COLLARLESS PRESS-FIT STANDARD OFFSET FEMORAL STEM	NO	METAL	HA COATED
Novation	1640208-1640218	HA COLLARLESS PRESS-FIT EXTENDED OFFSET FEMORAL STEM	NO	METAL	HA COATED
Novation	1640308-1640318	HA COLLARED PRESS-FIT STANDARD OFFSET FEMORAL STEM	NO	METAL	HA COATED
Novation	1641108-1641112	HA COLLARLESS PRESS-FIT STANDARD OFFSET FEMORAL STEM 12/14	NO	METAL	HA COATED
Novation	1641208-1641218	HA COLLARLESS PRESS-FIT EXTENDED OFFSET FEMORAL STEM	NO	METAL	HA COATED

Table 11: Revised Number of Novation Primary Total Conventional Hip Replacement by Catalogue Number Range

Femoral Stem Range	N Revised	N Total
1600009-1600018	0	10
1600109-1600118	2	3
1640108-1640118	55	912
1640208-1640218	31	353
1640308-1640318	10	101
1641108-1641112	2	32
1641208-1641218	0	12
TOTAL	100	1423

TABLE 12**Revision Rates of Novation Primary Total Conventional Hip Replacement by Component**

A prosthesis may be combined with multiple components. This analysis has been undertaken to determine if the revision rate varies according to the component with which it is combined.

Table 12: Revised Number of Novation Primary Total Conventional Hip Replacement by Acetabular Component

Acetabular Component	N Revised	N Total
Acetabular Shell (Global)	13	197
Adaptive	0	8
BSC-Cup	29	419
Cer-Met	1	2
Delta-TT	0	3
FMP	3	49
Fin II	1	16
Furlong	9	132
Logical G	17	231
Novae	5	99
Novae E	2	10
Novation	11	87
PINNACLE	1	41
Procotyl L	3	11
R3	0	1
Trident/Tritanium (Shell)	0	1
Trinity	4	115
ZCA	1	1
TOTAL	100	1423