

Dynasty Total Conventional Hip Investigation

Note: This analysis compares the Dynasty acetabular prosthesis with all other total conventional hip prostheses.

This prosthesis has been identified as having a significantly higher rate of revision. For a detailed explanation of the process used by the Registry that results in identification of prostheses that have a higher than anticipated rate of revision please refer to the Prostheses with Higher than Anticipated Rates of Revision chapter of the most recent AOANJRR Annual Report, <https://aoanjrr.sahmri.com/annual-reports-2026>.

Note: Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator. Procedures using prostheses with no recorded use in 2025 are excluded from the comparator.

TABLE 1
Revision Rate of Primary Total Conventional Hip Replacement

The revision rate of the Dynasty total conventional hip prosthesis is compared to all other total conventional hip prostheses.

Table 1: Revision Rates of Primary Total Conventional Hip Replacement

Component	N Revised	N Total	Obs. Years	Revisions/100 Obs. Yrs (95% CI)
Dynasty	123	2127	16103	0.76 (0.63, 0.91)
Other Total Conventional Hip	19701	568612	3623071	0.54 (0.54, 0.55)
TOTAL	19824	570739	3639174	0.54 (0.54, 0.55)

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 2

Yearly Cumulative Percent Revision of Primary Total Conventional Hip Replacement

The yearly cumulative percent revision of the Dynasty total conventional hip prosthesis is compared to all other total conventional hip prostheses.

Table 2: Yearly Cumulative Percent Revision (95% CI) of Primary Total Conventional Hip Replacement

CPR	1 Yr	2 Yrs	3 Yrs	4 Yrs	5 Yrs	6 Yrs
Dynasty	3.3 (2.7, 4.2)	3.9 (3.2, 4.8)	4.4 (3.6, 5.4)	4.6 (3.8, 5.6)	5.0 (4.1, 6.0)	5.3 (4.4, 6.4)
Other Total Conventional Hip	1.7 (1.7, 1.8)	2.1 (2.1, 2.2)	2.5 (2.4, 2.5)	2.8 (2.7, 2.8)	3.0 (3.0, 3.1)	3.3 (3.2, 3.3)

CPR	7 Yrs	8 Yrs	9 Yrs	10 Yrs	11 Yrs	12 Yrs
Dynasty	5.6 (4.7, 6.7)	5.9 (4.9, 7.0)	5.9 (4.9, 7.0)	6.2 (5.2, 7.5)	6.8 (5.6, 8.4)	6.8 (5.6, 8.4)
Other Total Conventional Hip	3.5 (3.5, 3.6)	3.8 (3.7, 3.9)	4.1 (4.0, 4.2)	4.4 (4.3, 4.5)	4.7 (4.6, 4.8)	5.1 (5.0, 5.2)

CPR	13 Yrs	14 Yrs	15 Yrs	16 Yrs	17 Yrs	18 Yrs
Dynasty	6.8 (5.6, 8.4)					
Other Total Conventional Hip	5.4 (5.3, 5.5)	5.8 (5.7, 5.9)	6.2 (6.1, 6.3)	6.6 (6.5, 6.7)	7.0 (6.8, 7.1)	7.3 (7.1, 7.5)

CPR	19 Yrs	20 Yrs	21 Yrs	22 Yrs	23 Yrs	24 Yrs	25 Yrs
Dynasty							
Other Total Conventional Hip	7.7 (7.5, 7.9)	8.0 (7.8, 8.2)	8.6 (8.3, 8.9)	9.0 (8.7, 9.4)	9.6 (9.1, 10.0)	9.6 (9.2, 10.1)	10.3 (9.4, 11.2)

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.


FIGURE 1

Yearly Cumulative Percent Revision of Primary Total Conventional Hip Replacement

The yearly cumulative percent revision of the Dynasty total conventional hip prosthesis is compared to all other total conventional hip prostheses. In addition, hazard ratios are reported.

Hazard ratios are reported for specific time periods during which the hazard ratio is constant. This is done to enable more specific and valid comparisons of the risk of revision over time. The pattern of variation in risk has important implications with respect to the underlying reasons for any difference.

Figure 1: Cumulative Percent Revision of Primary Total Conventional Hip Replacement

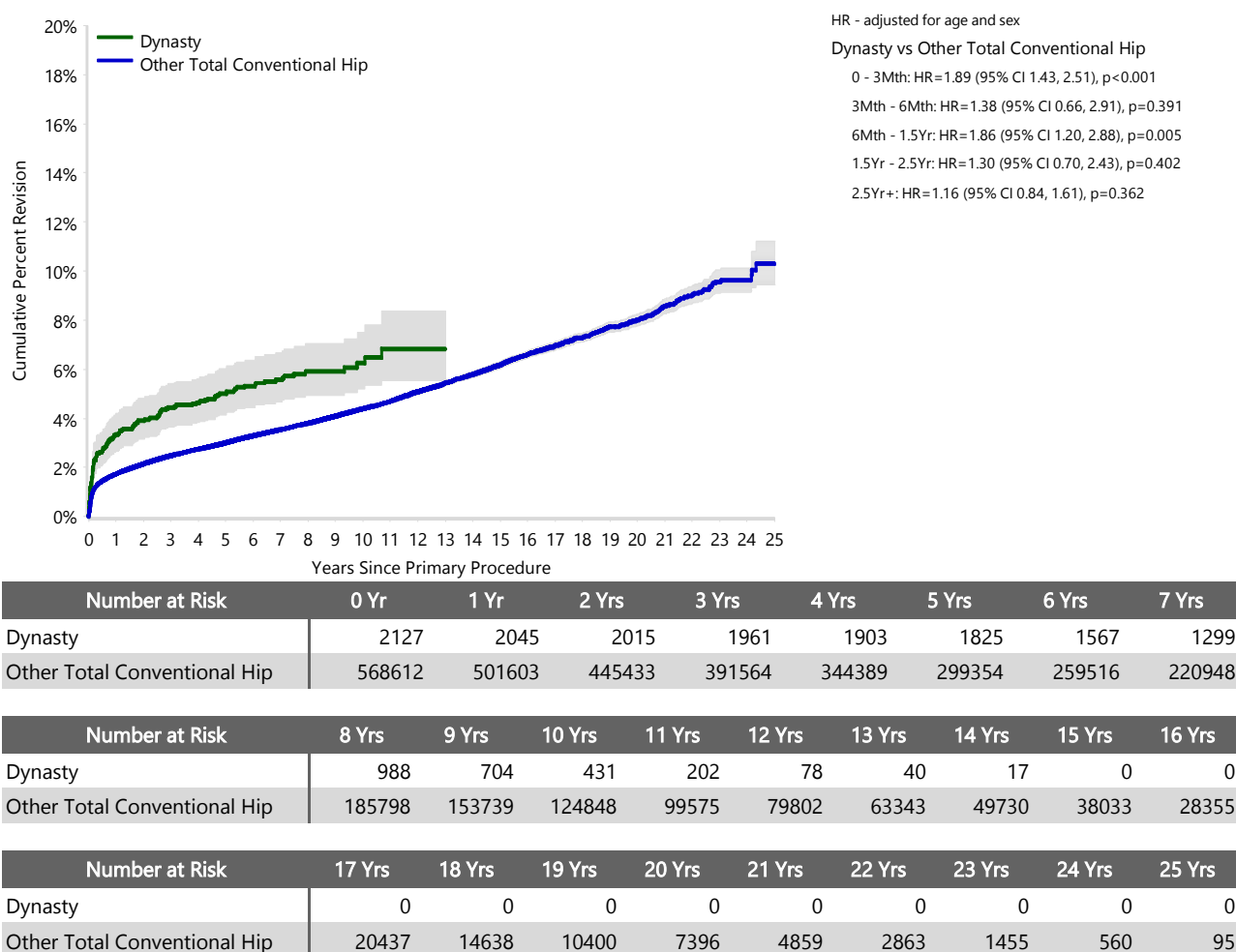

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 3**Primary Diagnosis for Revised Primary Total Conventional Hip Replacement**

This table identifies the diagnosis of the primary procedure which was subsequently revised. This information is provided as there is a variation on outcome depending on the primary diagnosis. It is therefore important when considering the reasons for a higher than anticipated rate of revision that there is identification of the primary diagnosis. This information should be compared to the primary diagnosis for the revisions of all other total conventional hip prostheses.

Table 3: Primary Diagnosis for Revised Primary Total Conventional Hip Replacement

Primary Diagnosis	Dynasty		Other Total Conventional Hip	
	Number	Percent	Number	Percent
Osteoarthritis	115	93.5	16296	82.7
Fractured Neck Of Femur	2	1.6	1521	7.7
Osteonecrosis	4	3.3	844	4.3
Developmental Dysplasia	2	1.6	321	1.6
Rheumatoid Arthritis			201	1.0
Failed Internal Fixation			165	0.8
Tumour			157	0.8
Other Inflammatory Arthritis			102	0.5
Fracture/Dislocation			60	0.3
Other			21	0.1
Arthrodesis Takedown			13	0.1
TOTAL	123	100.0	19701	100.0

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 4

Reasons for Revision

This is reported in two ways: a percentage of primary procedures revised and as a percentage of all revision procedures.

% Primaries Revised: This shows the proportional contribution of each revision diagnosis as a percentage of the total number of primary procedures. This percentage can be used to approximate the risk of being revised for that diagnosis. Differing percentages between groups, with the same distribution of follow up time, may identify problems of concern.

% Revisions: The number of revisions for each diagnosis is expressed as a percentage of the total number of revisions. This shows the distribution of reasons for revision within a group but cannot be used as a comparison between groups.

Table 4: Primary Total Conventional Hip Replacement - Reason for Revision (Follow-up Limited to 14.9 Years)

Revision Diagnosis	Number	Dynasty		Other Total Conventional Hip		
		% Primaries Revised	% Revisions	Number	% Primaries Revised	% Revisions
Infection	10	0.5	8.1	4959	0.9	25.8
Prosthesis Dislocation/Instability	29	1.4	23.6	4320	0.8	22.5
Fracture	39	1.8	31.7	4226	0.7	22.0
Loosening	28	1.3	22.8	3507	0.6	18.2
Pain	6	0.3	4.9	334	0.1	1.7
Leg Length Discrepancy	4	0.2	3.3	307	0.1	1.6
Malposition	2	0.1	1.6	280	0.0	1.5
Implant Breakage Stem				197	0.0	1.0
Lysis	1	0.0	0.8	171	0.0	0.9
Implant Breakage Acetabular Insert	1	0.0	0.8	124	0.0	0.6
Incorrect Sizing				96	0.0	0.5
Wear Acetabular Insert	2	0.1	1.6	85	0.0	0.4
Metal Related Pathology	1	0.0	0.8	74	0.0	0.4
Implant Breakage Acetabular				70	0.0	0.4
Tumour				46	0.0	0.2
Wear Head				45	0.0	0.2
Heterotopic Bone				26	0.0	0.1
Implant Breakage Head				23	0.0	0.1
Wear Acetabulum				6	0.0	0.0
Osteonecrosis				4	0.0	0.0
Synovitis				2	0.0	0.0
Progression Of Disease				1	0.0	0.0
Other				319	0.1	1.7
N Revision	123	5.8	100.0	19222	3.4	100.0
N Primary	2127			568612		

Note: This table is restricted to revisions within 14.9 years for all groups to allow a time-matched comparison of revisions.

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

FIGURE 2**Cumulative Incidence Revision Diagnosis of Primary Total Conventional Hip Replacement**

This figure details the cumulative incidence of the most common reasons for revision. The five most common reasons for revision are included as long as each of these reasons account for more than 10 procedures or at least 5% of all revisions for the Dynasty total conventional hip prosthesis. A comparative graph is provided of the cumulative incidence for the same reasons for revisions for all other total conventional hip prostheses.

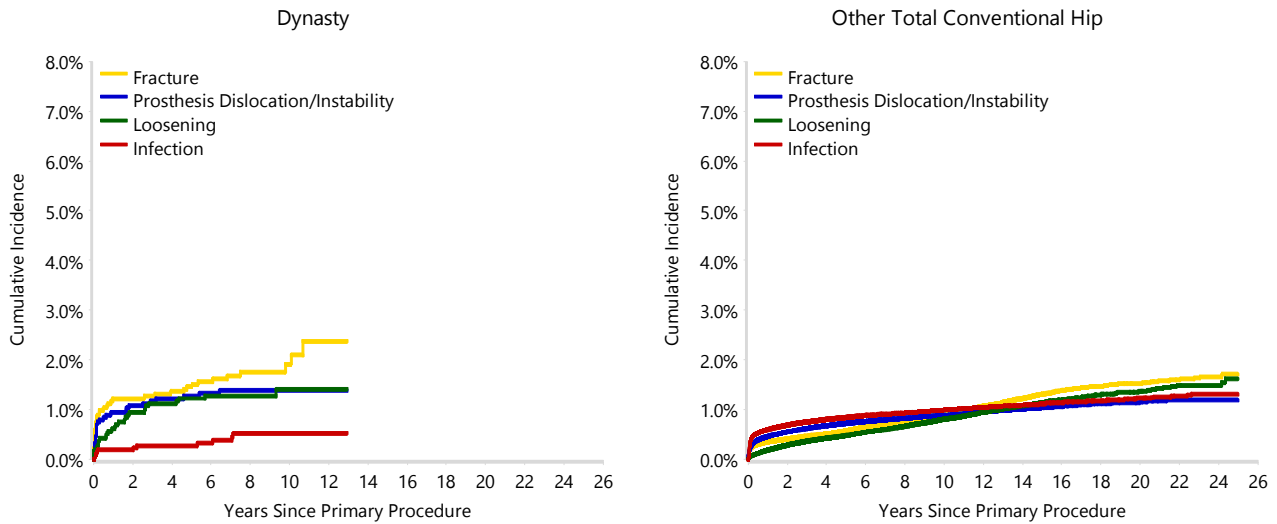
Figure 2: Cumulative Incidence Revision Diagnosis for Primary Total Conventional Hip Replacement

TABLE 5

Type of Revision Performed for Primary Total Conventional Hip Replacement

This analysis identifies the components used in the revision of the Dynasty total conventional hip prosthesis and compares it to the components used in the revision of all other total conventional hip prostheses.

The reason this analysis is undertaken is to identify whether there is one or more components which are being replaced that differ from the components replaced for revisions of all other total conventional hip prostheses i.e. is there a difference in the type of revision undertaken for the Dynasty total conventional hip prosthesis compared to all other total conventional hip prostheses.

Table 5: Primary Total Conventional Hip Replacement - Type of Revision (Follow-up Limited to 14.9 Years)

Type of Revision	Dynasty		Other Total Conventional Hip	
	Number	Percent	Number	Percent
Femoral Component	34	27.6	6480	33.7
Acetabular Component	48	39.0	3206	16.7
THR (Femoral/Acetabular)	15	12.2	2207	11.5
Cement Spacer	3	2.4	553	2.9
Removal of Prostheses	1	0.8	91	0.5
Reinsertion of Components	1	0.8	27	0.1
Total Femoral			12	0.1
Bipolar Head and Femoral			10	0.1
N Major	102	82.9	12586	65.5
Head/Insert	10	8.1	5285	27.5
Head Only	6	4.9	904	4.7
Minor Components	2	1.6	279	1.5
Insert Only	2	1.6	165	0.9
Bipolar Only			1	0.0
Cement Only			1	0.0
Head/Neck/Insert	1	0.8	1	0.0
N Minor	21	17.1	6636	34.5
TOTAL	123	100.0	19222	100.0

Note: This table is restricted to revisions within 14.9 years for all groups to allow a time-matched comparison of revisions.

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 6**Revision Rates of Dynasty Primary Total Conventional Hip Replacement by Fixation**

This analysis is provided as some prostheses have more than one fixation option. Additionally there are prostheses where an alternative to the recommended approach to fixation was used e.g. a cementless prosthesis that has been cemented or vice-versa.

Table 6: Revised Number of Dynasty Primary Total Conventional Hip Replacement by Fixation

Fixation	N Revised	N Total
Cementless	123	2113
Hybrid (Femur Cemented)	0	14
TOTAL	123	2127

TABLE 7**Revision Rates of Dynasty Primary Total Conventional Hip Replacement by Bearing Surface**

This analysis is provided as some prostheses are combined with a variety of bearing surfaces. All bearing surfaces used with this prosthesis are listed.

Table 7: Revised Number of Dynasty Primary Total Conventional Hip Replacement by Bearing Surface

Bearing Surface	N Revised	N Total
Ceramic/XLPE	83	1437
Metal/Metal	0	1
Metal/XLPE	39	686
Ceramicised Metal/XLPE	1	2
Unknown	0	1
TOTAL	123	2127

TABLE 8**Revision Rates of Dynasty Primary Total Conventional Hip Replacement by Approach**

This analysis is provided as some prostheses are used with a variety of surgical approaches. All surgical approaches used with this prosthesis are listed.

Table 8: Revised Number of Dynasty Primary Total Conventional Hip Replacement by Approach

Approach	N Revised	N Total
Anterior	40	368
Lateral	11	167
Posterior	45	1238
TOTAL	96	1773

Note: Excludes 354 procedures with no approach recorded

TABLE 9

Revision Rates of Primary Total Conventional Hip Replacement by State

This enables a state by state variation to be identified for the Dynasty total conventional hip prosthesis and provides the comparative data for each of the states for all other total conventional hip prostheses.

The purpose of this analysis is to determine if the higher than anticipated rate of revision has widespread distribution between states. If there is widespread distribution then the reason for the higher than anticipated rate of revision is unlikely to be surgeon specific. If the prosthesis has been used in only a small number of states it is not possible to distinguish if the higher than anticipated rate of revision is related to the prosthesis, surgeon, technique or patient.

Table 9: Revised Number of Primary Total Conventional Hip Replacement by State

Component	State	N Revised	N Total
Dynasty	NSW	103	1926
	VIC	0	1
	QLD	7	69
	WA	13	131
Other Total Conventional Hip	NSW	5507	167002
	VIC	4853	146558
	QLD	3968	104359
	WA	2510	64946
	SA	1887	52084
	TAS	439	18269
	ACT/NT	537	15394
TOTAL		19824	570739

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 10**Number of Revisions of Dynasty Primary Total Conventional Hip Replacement by Year of Implant**

This analysis details the number of prostheses reported each year to the Registry for the Dynasty total conventional hip prosthesis. It also provides the subsequent number of revisions of the primaries reported in that year.

Primary procedures performed in later years have had less follow up time therefore the number revised is expected to be less than the number revised in earlier years. For example, a primary procedure performed in 2025 has a maximum of one year to be revised, whereas a primary procedure performed in 2023 has a maximum of three years to be revised.

Table 10: Number of Revisions of Dynasty Primary Total Conventional Hip Replacement by Year of Implant

Year of Implant	Number Revised	Total Number
2011	3	40
2012	4	31
2013	3	49
2014	13	178
2015	19	298
2016	22	317
2017	17	306
2018	17	307
2019	14	272
2020	10	241
2021	1	39
2022	0	32
2023	0	17
TOTAL	123	2127

TABLE 11**Revision Rates of Dynasty Primary Total Conventional Hip Replacement by Catalogue Number Range**

Many prostheses have a number of catalogue ranges. The catalogue range is specific to particular design features; more than one catalogue range usually indicates a minor difference in design in a particular Dynasty prosthesis.

This analysis has been undertaken to determine if the revision rate varies according to the catalogue number range.

Model	Catalogue Range	Catalogue Description	Cement	Material
Acetabular				
Dynasty	DSBFGB46-DSBFGK76	DYNASTY BIOFOAM SHELL Ti6AL4V Ti	NO	METAL
Dynasty	DSPCGA42-DSPCGH68	DYNASTY PC SHELL Ti6AL4V	NO	METAL

Table 11: Revised Number of Dynasty Primary Total Conventional Hip Replacement by Catalogue Number Range

Acetabular Range	N Revised	N Total
DSBFGB46-DSBFGK76	99	1357
DSPCGA42-DSPCGH68	24	770
TOTAL	123	2127

TABLE 12

Revision Rates of Dynasty Primary Total Conventional Hip Replacement by Component

A prosthesis may be combined with multiple components. This analysis has been undertaken to determine if the revision rate varies according to the component with which it is combined.

Table 12: Revised Number of Dynasty Primary Total Conventional Hip Replacement by Femoral Stem Component

Femoral Stem Component	N Revised	N Total
CORAIL	0	2
CPCS	0	1
CPT	0	1
Exeter V40	0	9
Guardian	0	1
Profemur	0	1
Profemur E	0	11
Profemur Gladiator	6	105
Profemur L	110	1924
Profemur L (exch neck)	0	14
Profemur TL	7	56
Profemur XM	0	2
TOTAL	123	2127