

Genesis II Oxinium PS (cemented)/Genesis II (keel) Total Knee Investigation

Note: This analysis compares the Genesis II Oxinium PS (ctd)/Genesis II (keel) femoral/tibial combination with all other total knee prostheses.

This combination has been identified as having a significantly higher rate of revision. For a detailed explanation of the process used by the Registry that results in identification of prostheses that have a higher than anticipated rate of revision please refer to the Prostheses with Higher than Anticipated Rates of Revision chapter of the most recent AOANJRR Annual Report, <https://aoanjrr.sahmri.com/annual-reports-2026>.

Note: Procedures using prostheses with no recorded use in 2025 are excluded from the comparator.

TABLE 1

Revision Rate of Primary Total Knee Replacement

The revision rate of the Genesis II Oxinium PS (ctd)/Genesis II (keel) total knee combination is compared to all other total knee prostheses.

Table 1: Revision Rates of Primary Total Knee Replacement

Component	N Revised	N Total	Obs. Years	Revisions/100 Obs. Yrs (95% CI)
Genesis II Oxinium PS (ctd)/Genesis II (keel)	72	269	3363	2.14 (1.68, 2.70)
Other Total Knee	27212	791161	5142559	0.53 (0.52, 0.54)
TOTAL	27284	791430	5145922	0.53 (0.52, 0.54)

Note: Prostheses no longer used in 2025 are excluded from the comparator.

TABLE 2

Yearly Cumulative Percent Revision of Primary Total Knee Replacement

The yearly cumulative percent revision of the Genesis II Oxinium PS (ctd)/Genesis II (keel) total knee combination is compared to all other total knee prostheses.

Table 2: Yearly Cumulative Percent Revision (95% CI) of Primary Total Knee Replacement

CPR	1 Yr	2 Yrs	3 Yrs	4 Yrs	5 Yrs	6 Yrs
Genesis II Oxinium PS (ctd)/Genesis II (keel)	4.5 (2.6, 7.7)	11.2 (7.9, 15.6)	14.9 (11.1, 19.7)	17.1 (13.1, 22.2)	19.0 (14.8, 24.3)	19.8 (15.5, 25.1)
Other Total Knee	1.0 (1.0, 1.0)	1.8 (1.8, 1.8)	2.4 (2.3, 2.4)	2.8 (2.7, 2.8)	3.1 (3.0, 3.1)	3.4 (3.4, 3.4)

CPR	7 Yrs	8 Yrs	9 Yrs	10 Yrs	11 Yrs	12 Yrs
Genesis II Oxinium PS (ctd)/Genesis II (keel)	20.6 (16.2, 26.0)	22.2 (17.7, 27.7)	22.6 (18.0, 28.2)	22.6 (18.0, 28.2)	23.5 (18.9, 29.2)	25.4 (20.5, 31.2)
Other Total Knee	3.7 (3.7, 3.8)	4.0 (3.9, 4.1)	4.3 (4.2, 4.4)	4.6 (4.6, 4.7)	4.9 (4.9, 5.0)	5.2 (5.2, 5.3)

CPR	13 Yrs	14 Yrs	15 Yrs	16 Yrs	17 Yrs	18 Yrs
Genesis II Oxinium PS (ctd)/Genesis II (keel)	25.9 (21.0, 31.7)	25.9 (21.0, 31.7)	27.0 (21.9, 33.0)	27.0 (21.9, 33.0)	27.0 (21.9, 33.0)	28.4 (23.1, 34.6)
Other Total Knee	5.6 (5.5, 5.7)	5.9 (5.8, 6.0)	6.3 (6.2, 6.4)	6.7 (6.6, 6.8)	7.1 (7.0, 7.3)	7.5 (7.4, 7.7)

CPR	19 Yrs	20 Yrs	21 Yrs	22 Yrs	23 Yrs	24 Yrs	25 Yrs
Genesis II Oxinium PS (ctd)/Genesis II (keel)	29.1 (23.7, 35.4)						
Other Total Knee	7.8 (7.7, 8.0)	8.1 (7.9, 8.3)	8.5 (8.2, 8.7)	8.7 (8.4, 9.0)	8.8 (8.5, 9.1)	9.1 (8.7, 9.5)	9.5 (8.8, 10.2)

Note: Prostheses no longer used in 2025 are excluded from the comparator.


FIGURE 1

Yearly Cumulative Percent Revision of Primary Total Knee Replacement

The yearly cumulative percent revision of the Genesis II Oxinium PS (ctd)/Genesis II (keel) total knee combination is compared to all other total knee prostheses. In addition, hazard ratios are reported.

Hazard ratios are reported for specific time periods during which the hazard ratio is constant. This is done to enable more specific and valid comparisons of the risk of revision over time. The pattern of variation in risk has important implications with respect to the underlying reasons for any difference.

Figure 1: Cumulative Percent Revision of Primary Total Knee Replacement

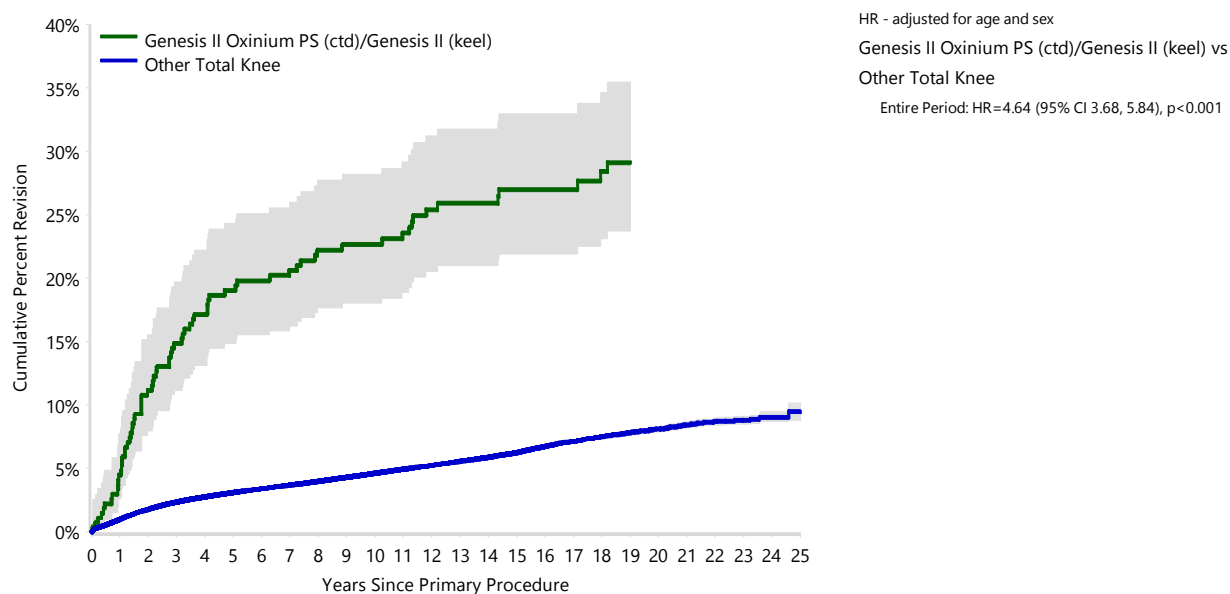

Number at Risk	0 Yr	1 Yr	2 Yrs	3 Yrs	4 Yrs	5 Yrs	6 Yrs	7 Yrs
Genesis II Oxinium PS (ctd)/Genesis II (keel)	269	257	239	227	219	212	205	199
Other Total Knee	791161	708442	628035	550608	487325	424656	369640	315800

Number at Risk	8 Yrs	9 Yrs	10 Yrs	11 Yrs	12 Yrs	13 Yrs	14 Yrs	15 Yrs	16 Yrs
Genesis II Oxinium PS (ctd)/Genesis II (keel)	191	185	173	167	152	146	140	130	123
Other Total Knee	266625	222135	182481	146813	117049	92041	70240	52396	38228

Number at Risk	17 Yrs	18 Yrs	19 Yrs	20 Yrs	21 Yrs	22 Yrs	23 Yrs	24 Yrs	25 Yrs
Genesis II Oxinium PS (ctd)/Genesis II (keel)	105	101	46	7	0	0	0	0	0
Other Total Knee	27523	19333	13570	9420	6164	3657	2096	972	244

Note: Prostheses no longer used in 2025 are excluded from the comparator.

TABLE 3**Primary Diagnosis for Revised Primary Total Knee Replacement**

This table identifies the diagnosis of the primary procedure which was subsequently revised. This information is provided as there is a variation on outcome depending on the primary diagnosis. It is therefore important when considering the reasons for a higher than anticipated rate of revision that there is identification of the primary diagnosis. This information should be compared to the primary diagnosis for the revisions of all other total knee prostheses.

Table 3: Primary Diagnosis for Revised Primary Total Knee Replacement

Primary Diagnosis	Genesis II Oxinium PS (ctd)/Genesis II (keel)		Other Total Knee	
	Number	Percent	Number	Percent
Osteoarthritis	71	98.6	26321	96.7
Rheumatoid Arthritis	1	1.4	329	1.2
Tumour			199	0.7
Other Inflammatory Arthritis			177	0.7
Osteonecrosis			104	0.4
Fracture			55	0.2
Other			26	0.1
Chondrocalcinosis			1	0.0
TOTAL	72	100.0	27212	100.0

Note: Prostheses no longer used in 2025 are excluded from the comparator.

TABLE 4

Reasons for Revision

This is reported in two ways: a percentage of primary procedures revised and as a percentage of all revision procedures.

% Primaries Revised: This shows the proportional contribution of each revision diagnosis as a percentage of the total number of primary procedures. This percentage can be used to approximate the risk of being revised for that diagnosis. Differing percentages between groups, with the same distribution of follow up time, may identify problems of concern.

% Revisions: The number of revisions for each diagnosis is expressed as a percentage of the total number of revisions. This shows the distribution of reasons for revision within a group but cannot be used as a comparison between groups.

Table 4: Primary Total Knee Replacement - Reason for Revision (Follow-up Limited to 20.4 Years)

Revision Diagnosis	Genesis II Oxinium PS (ctd)/Genesis II (keel)			Other Total Knee		
	Number	% Primaries Revised	% Revisions	Number	% Primaries Revised	% Revisions
Infection	12	4.5	16.7	7872	1.0	29.0
Loosening	50	18.6	69.4	5705	0.7	21.0
Instability				2871	0.4	10.6
Patella Erosion	2	0.7	2.8	2029	0.3	7.5
Pain	1	0.4	1.4	1865	0.2	6.9
Patellofemoral Pain				1556	0.2	5.7
Arthrofibrosis	1	0.4	1.4	1180	0.1	4.3
Fracture	1	0.4	1.4	1132	0.1	4.2
Malalignment				511	0.1	1.9
Wear Tibial Insert	2	0.7	2.8	370	0.0	1.4
Lysis				291	0.0	1.1
Implant Breakage Tibial Insert	2	0.7	2.8	250	0.0	0.9
Incorrect Sizing				218	0.0	0.8
Patella Maltracking				192	0.0	0.7
Bearing Dislocation				150	0.0	0.6
Implant Breakage Patella				147	0.0	0.5
Metal Related Pathology				114	0.0	0.4
Prosthesis Dislocation				78	0.0	0.3
Synovitis				64	0.0	0.2
Wear Patella				52	0.0	0.2
Osteonecrosis				51	0.0	0.2
Implant Breakage Femoral				49	0.0	0.2
Implant Breakage Tibial	1	0.4	1.4	38	0.0	0.1
Tumour				36	0.0	0.1
Heterotopic Bone				16	0.0	0.1
Progression Of Disease				8	0.0	0.0
Wear Tibial				6	0.0	0.0
Patella Dislocation				3	0.0	0.0
Incorrect Side				1	0.0	0.0
Wear Femoral				1	0.0	0.0
Other				314	0.0	1.2
N Revision	72	26.8	100.0	27170	3.4	100.0
N Primary	269			791161		

Note: This table is restricted to revisions within 20.4 years for all groups to allow a time-matched comparison of revisions.

Note: Prostheses no longer used in 2025 are excluded from the comparator.

FIGURE 2**Cumulative Incidence Revision Diagnosis of Primary Total Knee Replacement**

This figure details the cumulative incidence of the most common reasons for revision. The five most common reasons for revision are included as long as each of these reasons account for more than 10 procedures or at least 5% of all revisions for the Genesis II Oxinium PS (ctd)/Genesis II (keel) total knee combination. A comparative graph is provided of the cumulative incidence for the same reasons for revisions for all other total knee prostheses.

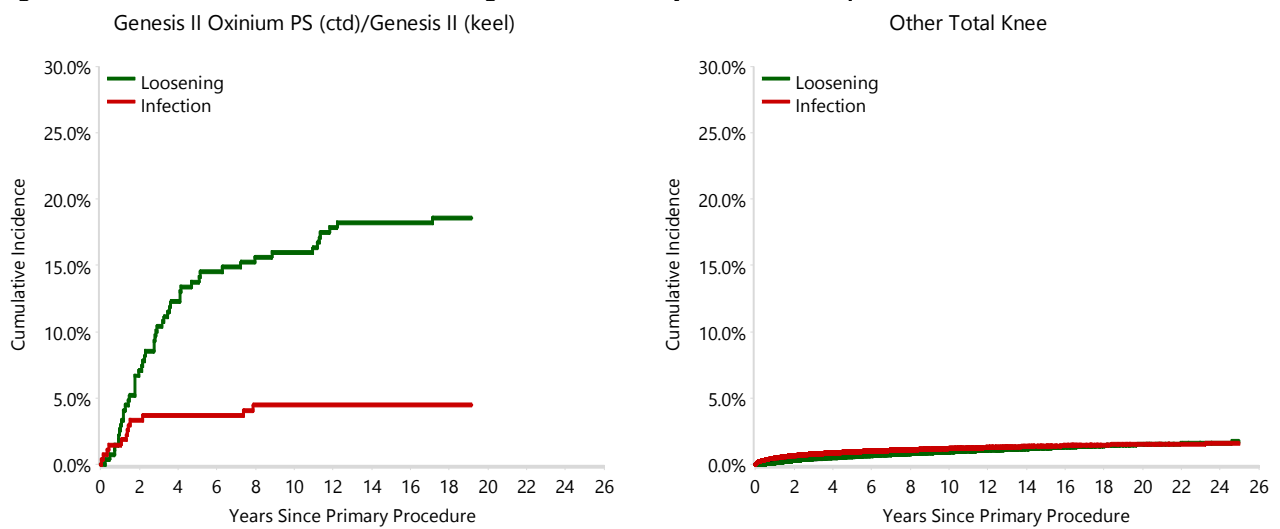
Figure 2: Cumulative Incidence Revision Diagnosis for Primary Total Knee Replacement

TABLE 5

Type of Revision Performed for Primary Total Knee Replacement

This analysis identifies the components used in the revision of the Genesis II Oxinium PS (ctd)/Genesis II (keel) total knee combination and compares it to the components used in the revision of all other total knee prostheses.

The reason this analysis is undertaken is to identify whether there is one or more components which are being replaced that differ from the components replaced for revisions of all other total knee prostheses i.e. is there a difference in the type of revision undertaken for the Genesis II Oxinium PS (ctd)/Genesis II (keel) total knee combination compared to all other total knee prostheses.

Table 5: Primary Total Knee Replacement - Type of Revision (Follow-up Limited to 20.4 Years)

Type of Revision	Genesis II Oxinium PS (ctd)/Genesis II (keel)		Other Total Knee	
	Number	Percent	Number	Percent
TKR (Tibial/Femoral)	11	15.3	6645	24.5
Tibial Component	41	56.9	1867	6.9
Femoral Component	2	2.8	1311	4.8
Cement Spacer	1	1.4	1136	4.2
Removal of Prostheses			134	0.5
Total Femoral			23	0.1
Reinsertion of Components			8	0.0
N Major	55	76.4	11124	40.9
Insert Only	10	13.9	8555	31.5
Patella Only	4	5.6	4488	16.5
Insert/Patella	3	4.2	2926	10.8
Minor Components			66	0.2
Cement Only			11	0.0
N Minor	17	23.6	16046	59.1
TOTAL	72	100.0	27170	100.0

Note: This table is restricted to revisions within 20.4 years for all groups to allow a time-matched comparison of revisions.

Note: Prostheses no longer used in 2025 are excluded from the comparator.

TABLE 6**Revision Rates of Genesis II Oxinium PS (ctd)/Genesis II (keel) Primary Total Knee Replacement by Fixation**

This analysis is provided as some prostheses have more than one fixation option. Additionally there are prostheses where an alternative to the recommended approach to fixation was used e.g. a cementless prosthesis that has been cemented or vice-versa.

Table 6: Revised Number of Genesis II Oxinium PS (ctd)/Genesis II (keel) Primary Total Knee Replacement by Fixation

Fixation	N Revised	N Total
Cemented	72	269
TOTAL	72	269

TABLE 7**Revision Rates of Genesis II Oxinium PS (ctd)/Genesis II (keel) Primary Total Knee Replacement by Bearing Surface**

This analysis is provided as some prostheses are combined with a variety of bearing surfaces. All bearing surfaces used with this combination are listed.

Table 7: Revised Number of Genesis II Oxinium PS (ctd)/Genesis II (keel) Primary Total Knee Replacement by Bearing Surface

Bearing Surface	N Revised	N Total
Non XLPE	71	268
Unknown	1	1
TOTAL	72	269

TABLE 8**Revision Rates of Genesis II Oxinium PS (ctd)/Genesis II (keel) Primary Total Knee Replacement by Bearing Mobility**

This analysis is provided as some prostheses are combined with a variety of bearing mobilities. All bearing mobilities used with this combination are listed.

Table 8: Revised Number of Genesis II Oxinium PS (ctd)/Genesis II (keel) Primary Total Knee Replacement by Bearing Mobility

Bearing Mobility	N Revised	N Total
Fixed	71	268
Unknown	1	1
TOTAL	72	269

TABLE 9**Revision Rates of Genesis II Oxinium PS (ctd)/Genesis II (keel) Primary Total Knee Replacement by Stability**

This analysis is provided as some prostheses are combined with a variety of stabilities. All stabilities used with this combination are listed.

Table 9: Revised Number of Genesis II Oxinium PS (ctd)/Genesis II (keel) Primary Total Knee Replacement by Stability

Stability	N Revised	N Total
Posterior Stabilised	71	268
Unknown	1	1
TOTAL	72	269

TABLE 10**Revision Rates of Primary Total Knee Replacement by State**

This enables a state by state variation to be identified for the Genesis II Oxinium PS (ctd)/Genesis II (keel) total knee combination and provides the comparative data for each of the states for all other total knee prostheses.

The purpose of this analysis is to determine if the higher than anticipated rate of revision has widespread distribution between states. If there is widespread distribution then the reason for the higher than anticipated rate of revision is unlikely to be surgeon specific. If the prosthesis has been used in only a small number of states it is not possible to distinguish if the higher than anticipated rate of revision is related to the prosthesis, surgeon, technique or patient.

Table 10: Revised Number of Primary Total Knee Replacement by State

Component	State	N Revised	N Total
Genesis II Oxinium PS (ctd)/Genesis II (keel)	NSW	50	124
	VIC	4	15
	QLD	7	80
	WA	11	50
Other Total Knee	NSW	7500	259373
	VIC	6016	163550
	QLD	5782	167164
	WA	3589	92791
	SA	3220	73114
	TAS	475	15075
	ACT/NT	630	20094
TOTAL		27284	791430

Note: Prostheses no longer used in 2025 are excluded from the comparator.

TABLE 11**Number of Revisions of Genesis II Oxinium PS (ctd)/Genesis II (keel) Primary Total Knee Replacement by Year of Implant**

This analysis details the number of prostheses reported each year to the Registry for the Genesis II Oxinium PS (ctd)/Genesis II (keel) total knee combination. It also provides the subsequent number of revisions of the primaries reported in that year.

Primary procedures performed in later years have had less follow up time therefore the number revised is expected to be less than the number revised in earlier years. For example, a primary procedure performed in 2025 has a maximum of one year to be revised, whereas a primary procedure performed in 2023 has a maximum of three years to be revised.

Table 11: Number of Revisions of Genesis II Oxinium PS (ctd)/Genesis II (keel) Primary Total Knee Replacement by Year of Implant

Year of Implant	Number Revised	Total Number
2005	4	19
2006	36	123
2007	32	127
TOTAL	72	269

TABLE 12

Revision Rates of Genesis II Oxinium PS (ctd)/Genesis II (keel) Primary Total Knee Replacement by Catalogue Number Range

Many prostheses have a number of catalogue ranges. The catalogue range is specific to particular design features; more than one catalogue range usually indicates a minor difference in design in a particular Genesis II Oxinium PS (ctd)/Genesis II (keel) prosthesis.

This analysis has been undertaken to determine if the revision rate varies according to the catalogue number range.

Model	Catalogue Range	Catalogue Description	Cement
Femoral			
Genesis II Oxinium PS	71421012-71421118	NONPOROUS PS OXINIUM FEMORAL COMPONENT	YES
Tibial			
Genesis II	71423102-71423118	NONPOROUS MINIMALLY INVASIVE TIBIAL BASEPLATE	YES

Table 12: Revised Number of Genesis II Oxinium PS (ctd)/Genesis II (keel) Primary Total Knee Replacement by Catalogue Number Range

Femoral Range	Tibial Range	N Revised	N Total
71421012-71421118	71423102-71423118	72	269
TOTAL		72	269