

Furlong Evolution Total Conventional Hip Investigation

Note: This analysis compares the Furlong Evolution femoral stem prosthesis with all other total conventional hip prostheses.

This prosthesis has been identified as having a significantly higher rate of revision. For a detailed explanation of the process used by the Registry that results in identification of prostheses that have a higher than anticipated rate of revision please refer to the Prostheses with Higher than Anticipated Rates of Revision chapter of the most recent AOANJRR Annual Report, <https://aoanjrr.sahmri.com/annual-reports-2026>.

Note: Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator. Procedures using prostheses with no recorded use in 2025 are excluded from the comparator.

TABLE 1

Revision Rate of Primary Total Conventional Hip Replacement

The revision rate of the Furlong Evolution total conventional hip prosthesis is compared to all other total conventional hip prostheses.

Table 1: Revision Rates of Primary Total Conventional Hip Replacement

Component	N Revised	N Total	Obs. Years	Revisions/100 Obs. Yrs (95% CI)
Furlong Evolution	43	521	3391	1.27 (0.92, 1.71)
Other Total Conventional Hip	19701	568612	3623071	0.54 (0.54, 0.55)
TOTAL	19744	569133	3626462	0.54 (0.54, 0.55)

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 2

Yearly Cumulative Percent Revision of Primary Total Conventional Hip Replacement

The yearly cumulative percent revision of the Furlong Evolution total conventional hip prosthesis is compared to all other total conventional hip prostheses.

Table 2: Yearly Cumulative Percent Revision (95% CI) of Primary Total Conventional Hip Replacement

CPR	1 Yr	2 Yrs	3 Yrs	4 Yrs	5 Yrs	6 Yrs
Furlong Evolution	5.2 (3.6, 7.5)	6.5 (4.7, 9.0)	7.3 (5.4, 9.9)	7.7 (5.7, 10.4)	8.4 (6.3, 11.1)	8.4 (6.3, 11.1)
Other Total Conventional Hip	1.7 (1.7, 1.8)	2.1 (2.1, 2.2)	2.5 (2.4, 2.5)	2.8 (2.7, 2.8)	3.0 (3.0, 3.1)	3.3 (3.2, 3.3)

CPR	7 Yrs	8 Yrs	9 Yrs	10 Yrs	11 Yrs	12 Yrs
Furlong Evolution	8.4 (6.3, 11.1)	8.4 (6.3, 11.1)	8.4 (6.3, 11.1)	8.4 (6.3, 11.1)	8.4 (6.3, 11.1)	
Other Total Conventional Hip	3.5 (3.5, 3.6)	3.8 (3.7, 3.9)	4.1 (4.0, 4.2)	4.4 (4.3, 4.5)	4.7 (4.6, 4.8)	5.1 (5.0, 5.2)

CPR	13 Yrs	14 Yrs	15 Yrs	16 Yrs	17 Yrs	18 Yrs
Furlong Evolution						
Other Total Conventional Hip	5.4 (5.3, 5.5)	5.8 (5.7, 5.9)	6.2 (6.1, 6.3)	6.6 (6.5, 6.7)	7.0 (6.8, 7.1)	7.3 (7.1, 7.5)

CPR	19 Yrs	20 Yrs	21 Yrs	22 Yrs	23 Yrs	24 Yrs	25 Yrs
Furlong Evolution							
Other Total Conventional Hip	7.7 (7.5, 7.9)	8.0 (7.8, 8.2)	8.6 (8.3, 8.9)	9.0 (8.7, 9.4)	9.6 (9.1, 10.0)	9.6 (9.2, 10.1)	10.3 (9.4, 11.2)

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.


FIGURE 1

Yearly Cumulative Percent Revision of Primary Total Conventional Hip Replacement

The yearly cumulative percent revision of the Furlong Evolution total conventional hip prosthesis is compared to all other total conventional hip prostheses. In addition, hazard ratios are reported.

Hazard ratios are reported for specific time periods during which the hazard ratio is constant. This is done to enable more specific and valid comparisons of the risk of revision over time. The pattern of variation in risk has important implications with respect to the underlying reasons for any difference.

Figure 1: Cumulative Percent Revision of Primary Total Conventional Hip Replacement

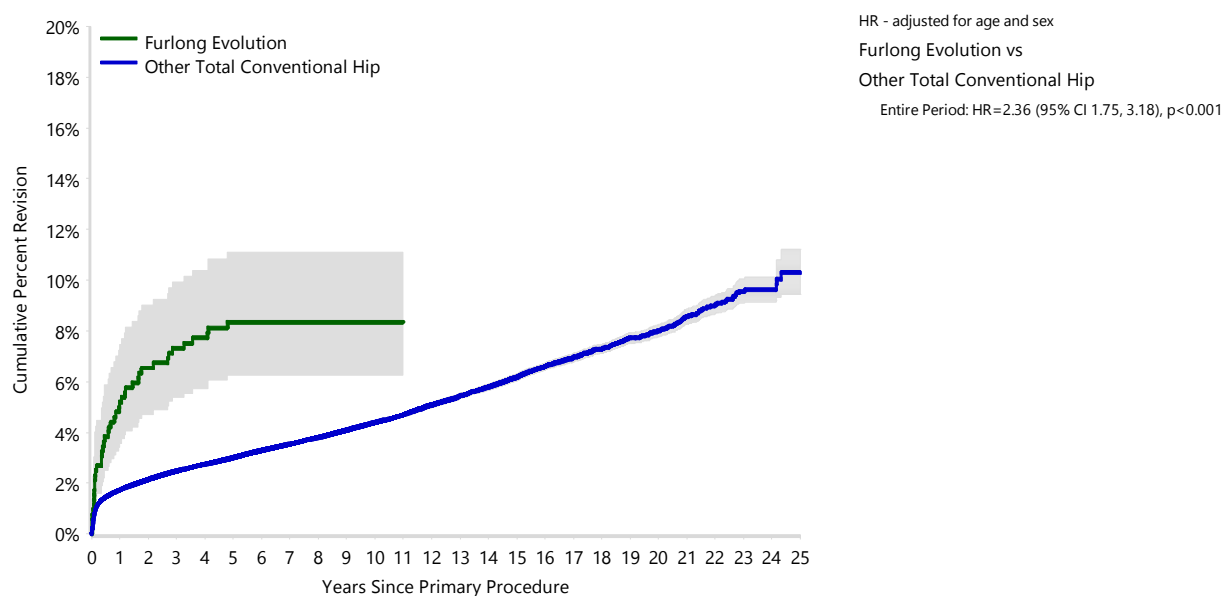

Number at Risk	0 Yr	1 Yr	2 Yrs	3 Yrs	4 Yrs	5 Yrs	6 Yrs	7 Yrs
Furlong Evolution	521	492	482	470	463	385	308	212
Other Total Conventional Hip	568612	501603	445433	391564	344389	299354	259516	220948

Number at Risk	8 Yrs	9 Yrs	10 Yrs	11 Yrs	12 Yrs	13 Yrs	14 Yrs	15 Yrs	16 Yrs
Furlong Evolution	123	78	67	43	23	0	0	0	0
Other Total Conventional Hip	185798	153739	124848	99575	79802	63343	49730	38033	28355

Number at Risk	17 Yrs	18 Yrs	19 Yrs	20 Yrs	21 Yrs	22 Yrs	23 Yrs	24 Yrs	25 Yrs
Furlong Evolution	0	0	0	0	0	0	0	0	0
Other Total Conventional Hip	20437	14638	10400	7396	4859	2863	1455	560	95

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 3**Primary Diagnosis for Revised Primary Total Conventional Hip Replacement**

This table identifies the diagnosis of the primary procedure which was subsequently revised. This information is provided as there is a variation on outcome depending on the primary diagnosis. It is therefore important when considering the reasons for a higher than anticipated rate of revision that there is identification of the primary diagnosis. This information should be compared to the primary diagnosis for the revisions of all other total conventional hip prostheses.

Table 3: Primary Diagnosis for Revised Primary Total Conventional Hip Replacement

Primary Diagnosis	Furlong Evolution		Other Total Conventional Hip	
	Number	Percent	Number	Percent
Osteoarthritis	40	93.0	16296	82.7
Fractured Neck Of Femur	1	2.3	1521	7.7
Osteonecrosis	2	4.7	844	4.3
Developmental Dysplasia			321	1.6
Rheumatoid Arthritis			201	1.0
Failed Internal Fixation			165	0.8
Tumour			157	0.8
Other Inflammatory Arthritis			102	0.5
Fracture/Dislocation			60	0.3
Other			21	0.1
Arthrodesis Takedown			13	0.1
TOTAL	43	100.0	19701	100.0

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 4

Reasons for Revision

This is reported in two ways: a percentage of primary procedures revised and as a percentage of all revision procedures.

% Primaries Revised: This shows the proportional contribution of each revision diagnosis as a percentage of the total number of primary procedures. This percentage can be used to approximate the risk of being revised for that diagnosis. Differing percentages between groups, with the same distribution of follow up time, may identify problems of concern.

% Revisions: The number of revisions for each diagnosis is expressed as a percentage of the total number of revisions. This shows the distribution of reasons for revision within a group but cannot be used as a comparison between groups.

Table 4: Primary Total Conventional Hip Replacement - Reason for Revision (Follow-up Limited to 12.6 Years)

Revision Diagnosis	Furlong Evolution			Other Total Conventional Hip		
	Number	% Primaries Revised	% Revisions	Number	% Primaries Revised	% Revisions
Infection	12	2.3	27.9	4912	0.9	26.2
Prosthesis Dislocation/Instability	3	0.6	7.0	4269	0.8	22.8
Fracture	6	1.2	14.0	4059	0.7	21.6
Loosening	14	2.7	32.6	3377	0.6	18.0
Pain				327	0.1	1.7
Leg Length Discrepancy	3	0.6	7.0	306	0.1	1.6
Malposition	4	0.8	9.3	280	0.0	1.5
Implant Breakage Stem				187	0.0	1.0
Lysis				152	0.0	0.8
Implant Breakage Acetabular Insert				119	0.0	0.6
Incorrect Sizing	1	0.2	2.3	96	0.0	0.5
Implant Breakage Acetabular				69	0.0	0.4
Metal Related Pathology				68	0.0	0.4
Wear Acetabular Insert				65	0.0	0.3
Tumour				45	0.0	0.2
Wear Head				42	0.0	0.2
Heterotopic Bone				26	0.0	0.1
Implant Breakage Head				22	0.0	0.1
Wear Acetabulum				5	0.0	0.0
Osteonecrosis				4	0.0	0.0
Synovitis				2	0.0	0.0
Progression Of Disease				1	0.0	0.0
Other				317	0.1	1.7
N Revision	43	8.3	100.0	18750	3.3	100.0
N Primary	521			568612		

Note: This table is restricted to revisions within 12.6 years for all groups to allow a time-matched comparison of revisions.

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

FIGURE 2**Cumulative Incidence Revision Diagnosis of Primary Total Conventional Hip Replacement**

This figure details the cumulative incidence of the most common reasons for revision. The five most common reasons for revision are included as long as each of these reasons account for more than 10 procedures or at least 5% of all revisions for the Furlong Evolution total conventional hip prosthesis. A comparative graph is provided of the cumulative incidence for the same reasons for revisions for all other total conventional hip prostheses.

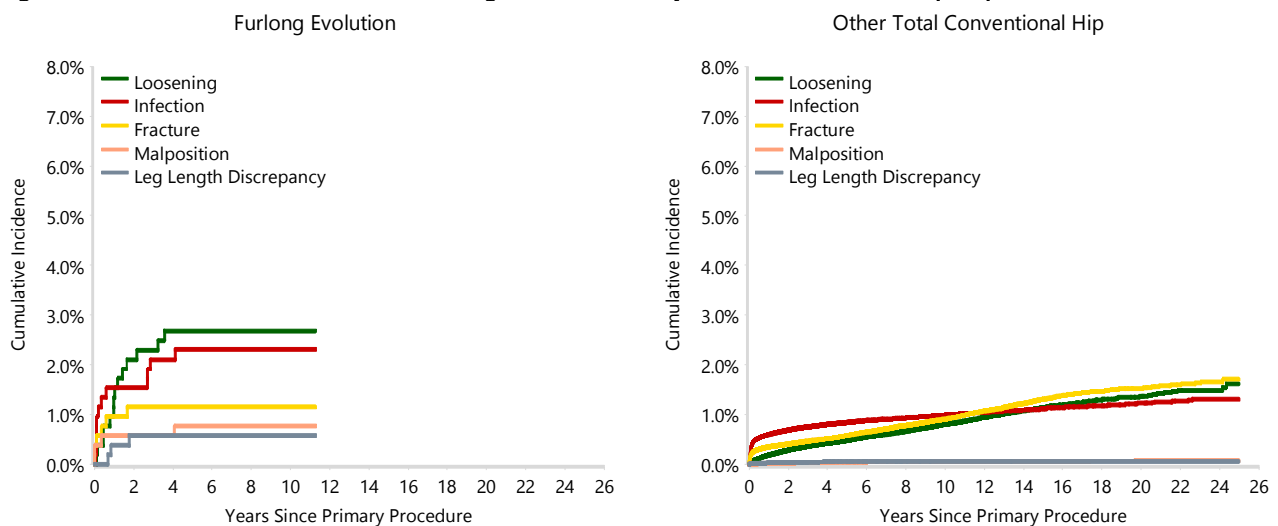
Figure 2: Cumulative Incidence Revision Diagnosis for Primary Total Conventional Hip Replacement

TABLE 5

Type of Revision Performed for Primary Total Conventional Hip Replacement

This analysis identifies the components used in the revision of the Furlong Evolution total conventional hip prosthesis and compares it to the components used in the revision of all other total conventional hip prostheses.

The reason this analysis is undertaken is to identify whether there is one or more components which are being replaced that differ from the components replaced for revisions of all other total conventional hip prostheses i.e. is there a difference in the type of revision undertaken for the Furlong Evolution total conventional hip prosthesis compared to all other total conventional hip prostheses.

Table 5: Primary Total Conventional Hip Replacement - Type of Revision (Follow-up Limited to 12.6 Years)

Type of Revision	Furlong Evolution		Other Total Conventional Hip	
	Number	Percent	Number	Percent
Femoral Component	17	39.5	6261	33.4
Acetabular Component	8	18.6	3118	16.6
THR (Femoral/Acetabular)	6	14.0	2135	11.4
Cement Spacer	3	7.0	547	2.9
Removal of Prostheses	1	2.3	90	0.5
Reinsertion of Components			27	0.1
Total Femoral			11	0.1
Bipolar Head and Femoral			10	0.1
N Major	35	81.4	12199	65.1
Head/Insert	4	9.3	5210	27.8
Head Only	3	7.0	899	4.8
Minor Components	1	2.3	274	1.5
Insert Only			165	0.9
Bipolar Only			1	0.0
Cement Only			1	0.0
Head/Neck/Insert			1	0.0
N Minor	8	18.6	6551	34.9
TOTAL	43	100.0	18750	100.0

Note: This table is restricted to revisions within 12.6 years for all groups to allow a time-matched comparison of revisions.

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 6**Revision Rates of Furlong Evolution Primary Total Conventional Hip Replacement by Fixation**

This analysis is provided as some prostheses have more than one fixation option. Additionally there are prostheses where an alternative to the recommended approach to fixation was used e.g. a cementless prosthesis that has been cemented or vice-versa.

Table 6: Revised Number of Furlong Evolution Primary Total Conventional Hip Replacement by Fixation

Fixation	N Revised	N Total
Cementless	43	521
TOTAL	43	521

TABLE 7**Revision Rates of Furlong Evolution Primary Total Conventional Hip Replacement by Bearing Surface**

This analysis is provided as some prostheses are combined with a variety of bearing surfaces. All bearing surfaces used with this prosthesis are listed.

Table 7: Revised Number of Furlong Evolution Primary Total Conventional Hip Replacement by Bearing Surface

Bearing Surface	N Revised	N Total
Ceramic/Ceramic	23	253
Ceramic/Non XLPE	14	211
Ceramic/XLPE	1	17
Ceramic/XLPE + Antioxidant	0	4
Metal/Non XLPE	4	35
Metal/XLPE	1	1
TOTAL	43	521

TABLE 8**Revision Rates of Furlong Evolution Primary Total Conventional Hip Replacement by Approach**

This analysis is provided as some prostheses are used with a variety of surgical approaches. All surgical approaches used with this prosthesis are listed.

Table 8: Revised Number of Furlong Evolution Primary Total Conventional Hip Replacement by Approach

Approach	N Revised	N Total
Anterior	23	271
Lateral	5	27
Posterior	9	162
TOTAL	37	460

Note: Excludes 61 procedures with no approach recorded

TABLE 9

Revision Rates of Primary Total Conventional Hip Replacement by State

This enables a state by state variation to be identified for the Furlong Evolution total conventional hip prosthesis and provides the comparative data for each of the states for all other total conventional hip prostheses.

The purpose of this analysis is to determine if the higher than anticipated rate of revision has widespread distribution between states. If there is widespread distribution then the reason for the higher than anticipated rate of revision is unlikely to be surgeon specific. If the prosthesis has been used in only a small number of states it is not possible to distinguish if the higher than anticipated rate of revision is related to the prosthesis, surgeon, technique or patient.

Table 9: Revised Number of Primary Total Conventional Hip Replacement by State

Component	State	N Revised	N Total
Furlong Evolution	NSW	4	35
	VIC	13	201
	QLD	5	69
	WA	21	211
	SA	0	2
	ACT/NT	0	3
Other Total Conventional Hip	NSW	5507	167002
	VIC	4853	146558
	QLD	3968	104359
	WA	2510	64946
	SA	1887	52084
	TAS	439	18269
	ACT/NT	537	15394
TOTAL		19744	569133

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 10**Number of Revisions of Furlong Evolution Primary Total Conventional Hip Replacement by Year of Implant**

This analysis details the number of prostheses reported each year to the Registry for the Furlong Evolution total conventional hip prosthesis. It also provides the subsequent number of revisions of the primaries reported in that year.

Primary procedures performed in later years have had less follow up time therefore the number revised is expected to be less than the number revised in earlier years. For example, a primary procedure performed in 2025 has a maximum of one year to be revised, whereas a primary procedure performed in 2023 has a maximum of three years to be revised.

Table 10: Number of Revisions of Furlong Evolution Primary Total Conventional Hip Replacement by Year of Implant

Year of Implant	Number Revised	Total Number
2013	2	29
2014	4	25
2015	2	32
2016	0	11
2017	3	54
2018	11	102
2019	11	106
2020	7	83
2021	1	75
2022	2	4
TOTAL	43	521

TABLE 11

Revision Rates of Furlong Evolution Primary Total Conventional Hip Replacement by Catalogue Number Range

Many prostheses have a number of catalogue ranges. The catalogue range is specific to particular design features; more than one catalogue range usually indicates a minor difference in design in a particular Furlong Evolution prosthesis.

This analysis has been undertaken to determine if the revision rate varies according to the catalogue number range.

Model	Catalogue Range	Catalogue Description	Cement	Material
Femoral Stem				
Furlong Evolution	4260106-4260117	FURLONG EVOLUTION COLLARED 126° NECK ANGLE STD OFFSET	NO	METAL
Furlong Evolution	4260210-4260212	FURLONG EVOLUTION COLLARLESS 126° NECK ANGLE STD OFFSET Ti6AL4V	NO	METAL
Furlong Evolution	4265106-4265117	FURLONG EVOLUTION COLLARED 126° NECK ANGLE HIGH OFFSET Ti6AL4V	NO	METAL
Furlong Evolution	4265111-4265211	FURLONG EVOLUTION COLLARLESS 126° NECK ANGLE HIGH OFFSET Ti6AL4V	NO	METAL
Furlong Evolution	4330106-4330117	FURLONG EVOLUTION COLLARED 133° NECK ANGLE STD OFFSET Ti6AL4V	NO	METAL
Furlong Evolution	4330206-4330214	FURLONG EVOLUTION COLLARLESS 133° NECK ANGLE STD OFFSET Ti6AL4V	NO	METAL
Furlong Evolution	4335106-4335117	FURLONG EVOLUTION COLLARED 133° NECK ANGLE HIGH OFFSET Ti6AL4V	NO	METAL
Furlong Evolution	4335208-4335217	FURLONG EVOLUTION COLLARLESS 133° NECK ANGLE HIGH OFFSET Ti6AL4V	NO	METAL

Table 11: Revised Number of Furlong Evolution Primary Total Conventional Hip Replacement by Catalogue Number Range

Femoral Stem Range	N Revised	N Total
4260106-4260117	8	39
4260210-4260212	0	2
4265106-4265117	5	45
4265111-4265211	1	15
4330106-4330117	18	244
4330206-4330214	4	39
4335106-4335117	4	115
4335208-4335217	3	22
TOTAL	43	521

TABLE 12**Revision Rates of Furlong Evolution Primary Total Conventional Hip Replacement by Component**

A prosthesis may be combined with multiple components. This analysis has been undertaken to determine if the revision rate varies according to the component with which it is combined.

Table 12: Revised Number of Furlong Evolution Primary Total Conventional Hip Replacement by Acetabular Component

Acetabular Component	N Revised	N Total
Acetabular Shell (Global)	1	6
Adaptive	0	9
Delta-One-TT	0	1
Delta-TT	0	3
Furlong	37	437
Logical G	0	17
Mpact	4	41
PINNACLE	0	2
R3	1	1
Trident/Tritanium (Shell)	0	1
Versafitcup CC	0	3
TOTAL	43	521