

Adapter (cementless) Total Conventional Hip Investigation

Note: This analysis compares the Adapter (cless) femoral stem prosthesis with all other total conventional hip prostheses.

This prosthesis has been identified as having a significantly higher rate of revision. For a detailed explanation of the process used by the Registry that results in identification of prostheses that have a higher than anticipated rate of revision please refer to the Prostheses with Higher than Anticipated Rates of Revision chapter of the most recent AOANJRR Annual Report, <https://aoanjrr.sahmri.com/annual-reports-2026>.

Note: Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator. Procedures using prostheses with no recorded use in 2025 are excluded from the comparator.

TABLE 1

Revision Rate of Primary Total Conventional Hip Replacement

The revision rate of the Adapter (cless) total conventional hip prosthesis is compared to all other total conventional hip prostheses.

Table 1: Revision Rates of Primary Total Conventional Hip Replacement

Component	N Revised	N Total	Obs. Years	Revisions/100 Obs. Yrs (95% CI)
Adapter (cless)	162	744	8077	2.01 (1.71, 2.34)
Other Total Conventional Hip	19701	568612	3623071	0.54 (0.54, 0.55)
TOTAL	19863	569356	3631148	0.55 (0.54, 0.55)

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 2

Yearly Cumulative Percent Revision of Primary Total Conventional Hip Replacement

The yearly cumulative percent revision of the Adapter (cless) total conventional hip prosthesis is compared to all other total conventional hip prostheses.

Table 2: Yearly Cumulative Percent Revision (95% CI) of Primary Total Conventional Hip Replacement

CPR	1 Yr	2 Yrs	3 Yrs	4 Yrs	5 Yrs	6 Yrs
Adapter (cless)	3.2 (2.2, 4.8)	4.9 (3.6, 6.7)	6.9 (5.2, 8.9)	9.5 (7.5, 11.8)	11.7 (9.5, 14.3)	13.8 (11.5, 16.6)
Other Total Conventional Hip	1.7 (1.7, 1.8)	2.1 (2.1, 2.2)	2.5 (2.4, 2.5)	2.8 (2.7, 2.8)	3.0 (3.0, 3.1)	3.3 (3.2, 3.3)

CPR	7 Yrs	8 Yrs	9 Yrs	10 Yrs	11 Yrs	12 Yrs
Adapter (cless)	16.1 (13.5, 19.0)	18.1 (15.4, 21.2)	18.9 (16.2, 22.1)	20.1 (17.2, 23.3)	21.7 (18.7, 25.1)	23.2 (20.0, 26.7)
Other Total Conventional Hip	3.5 (3.5, 3.6)	3.8 (3.7, 3.9)	4.1 (4.0, 4.2)	4.4 (4.3, 4.5)	4.7 (4.6, 4.8)	5.1 (5.0, 5.2)

CPR	13 Yrs	14 Yrs	15 Yrs	16 Yrs	17 Yrs	18 Yrs
Adapter (cless)	24.5 (21.3, 28.2)	25.2 (21.9, 29.0)	25.8 (22.4, 29.5)	25.8 (22.4, 29.5)	26.2 (22.8, 30.0)	26.2 (22.8, 30.0)
Other Total Conventional Hip	5.4 (5.3, 5.5)	5.8 (5.7, 5.9)	6.2 (6.1, 6.3)	6.6 (6.5, 6.7)	7.0 (6.8, 7.1)	7.3 (7.1, 7.5)

CPR	19 Yrs	20 Yrs	21 Yrs	22 Yrs	23 Yrs	24 Yrs	25 Yrs
Adapter (cless)	27.0 (23.3, 31.2)						
Other Total Conventional Hip	7.7 (7.5, 7.9)	8.0 (7.8, 8.2)	8.6 (8.3, 8.9)	9.0 (8.7, 9.4)	9.6 (9.1, 10.0)	9.6 (9.2, 10.1)	10.3 (9.4, 11.2)

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.


FIGURE 1

Yearly Cumulative Percent Revision of Primary Total Conventional Hip Replacement

The yearly cumulative percent revision of the Adapter (cless) total conventional hip prosthesis is compared to all other total conventional hip prostheses. In addition, hazard ratios are reported.

Hazard ratios are reported for specific time periods during which the hazard ratio is constant. This is done to enable more specific and valid comparisons of the risk of revision over time. The pattern of variation in risk has important implications with respect to the underlying reasons for any difference.

Figure 1: Cumulative Percent Revision of Primary Total Conventional Hip Replacement

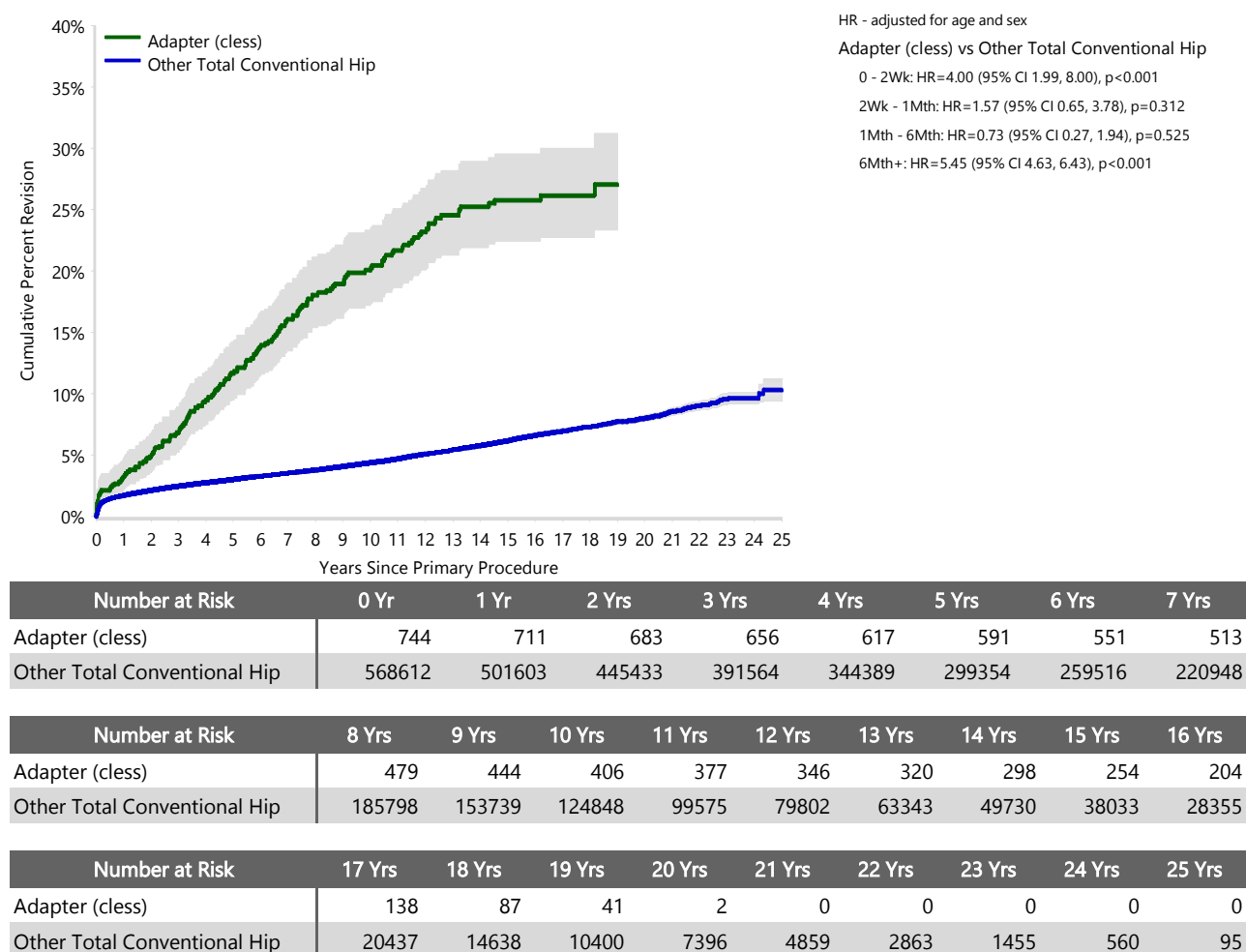

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 3**Primary Diagnosis for Revised Primary Total Conventional Hip Replacement**

This table identifies the diagnosis of the primary procedure which was subsequently revised. This information is provided as there is a variation on outcome depending on the primary diagnosis. It is therefore important when considering the reasons for a higher than anticipated rate of revision that there is identification of the primary diagnosis. This information should be compared to the primary diagnosis for the revisions of all other total conventional hip prostheses.

Table 3: Primary Diagnosis for Revised Primary Total Conventional Hip Replacement

Primary Diagnosis	Adapter (class)		Other Total Conventional Hip	
	Number	Percent	Number	Percent
Osteoarthritis	143	88.3	16296	82.7
Fractured Neck Of Femur	6	3.7	1521	7.7
Osteonecrosis	7	4.3	844	4.3
Developmental Dysplasia	1	0.6	321	1.6
Rheumatoid Arthritis	2	1.2	201	1.0
Failed Internal Fixation	1	0.6	165	0.8
Tumour			157	0.8
Other Inflammatory Arthritis	2	1.2	102	0.5
Fracture/Dislocation			60	0.3
Other			21	0.1
Arthrodesis Takedown			13	0.1
TOTAL	162	100.0	19701	100.0

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 4

Reasons for Revision

This is reported in two ways: a percentage of primary procedures revised and as a percentage of all revision procedures.

% Primaries Revised: This shows the proportional contribution of each revision diagnosis as a percentage of the total number of primary procedures. This percentage can be used to approximate the risk of being revised for that diagnosis. Differing percentages between groups, with the same distribution of follow up time, may identify problems of concern.

% Revisions: The number of revisions for each diagnosis is expressed as a percentage of the total number of revisions. This shows the distribution of reasons for revision within a group but cannot be used as a comparison between groups.

Table 4: Primary Total Conventional Hip Replacement - Reason for Revision (Follow-up Limited to 20.4 Years)

Revision Diagnosis	Number	Adapter (class)		Other Total Conventional Hip		
		% Primaries Revised	% Revisions	Number	% Primaries Revised	% Revisions
Infection	27	3.6	16.7	5010	0.9	25.5
Prosthesis Dislocation/Instability	23	3.1	14.2	4379	0.8	22.3
Fracture	5	0.7	3.1	4338	0.8	22.1
Loosening	33	4.4	20.4	3613	0.6	18.4
Pain	5	0.7	3.1	339	0.1	1.7
Leg Length Discrepancy	1	0.1	0.6	307	0.1	1.6
Malposition	3	0.4	1.9	282	0.0	1.4
Implant Breakage Stem	2	0.3	1.2	208	0.0	1.1
Lysis	16	2.2	9.9	196	0.0	1.0
Implant Breakage Acetabular Insert	1	0.1	0.6	126	0.0	0.6
Wear Acetabular Insert	2	0.3	1.2	106	0.0	0.5
Incorrect Sizing	1	0.1	0.6	96	0.0	0.5
Metal Related Pathology	37	5.0	22.8	90	0.0	0.5
Implant Breakage Acetabular	3	0.4	1.9	73	0.0	0.4
Wear Head				48	0.0	0.2
Tumour				46	0.0	0.2
Heterotopic Bone				26	0.0	0.1
Implant Breakage Head				26	0.0	0.1
Wear Acetabulum	1	0.1	0.6	6	0.0	0.0
Osteonecrosis				4	0.0	0.0
Synovitis				2	0.0	0.0
Progression Of Disease				1	0.0	0.0
Other	2	0.3	1.2	321	0.1	1.6
N Revision	162	21.8	100.0	19643	3.5	100.0
N Primary	744			568612		

Note: This table is restricted to revisions within 20.4 years for all groups to allow a time-matched comparison of revisions.

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

FIGURE 2**Cumulative Incidence Revision Diagnosis of Primary Total Conventional Hip Replacement**

This figure details the cumulative incidence of the most common reasons for revision. The five most common reasons for revision are included as long as each of these reasons account for more than 10 procedures or at least 5% of all revisions for the Adapter (cless) total conventional hip prosthesis. A comparative graph is provided of the cumulative incidence for the same reasons for revisions for all other total conventional hip prostheses.

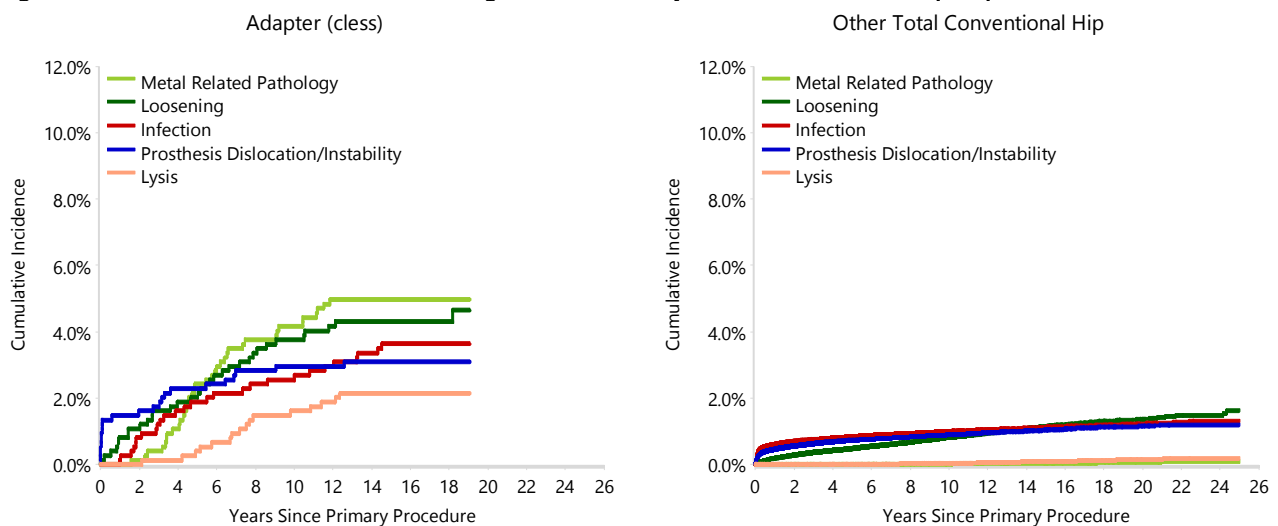
Figure 2: Cumulative Incidence Revision Diagnosis for Primary Total Conventional Hip Replacement

TABLE 5

Type of Revision Performed for Primary Total Conventional Hip Replacement

This analysis identifies the components used in the revision of the Adapter (cless) total conventional hip prosthesis and compares it to the components used in the revision of all other total conventional hip prostheses.

The reason this analysis is undertaken is to identify whether there is one or more components which are being replaced that differ from the components replaced for revisions of all other total conventional hip prostheses i.e. is there a difference in the type of revision undertaken for the Adapter (cless) total conventional hip prosthesis compared to all other total conventional hip prostheses.

Table 5: Primary Total Conventional Hip Replacement - Type of Revision (Follow-up Limited to 20.4 Years)

Type of Revision	Adapter (cless)		Other Total Conventional Hip	
	Number	Percent	Number	Percent
Femoral Component	15	9.3	6633	33.8
Acetabular Component	52	32.1	3297	16.8
THR (Femoral/Acetabular)	33	20.4	2288	11.6
Cement Spacer	9	5.6	556	2.8
Removal of Prostheses			92	0.5
Reinsertion of Components			27	0.1
Total Femoral			12	0.1
Bipolar Head and Femoral			10	0.1
N Major	109	67.3	12915	65.7
Head/Insert	4	2.5	5360	27.3
Head Only	3	1.9	912	4.6
Minor Components	3	1.9	285	1.5
Insert Only	1	0.6	168	0.9
Head/Neck/Insert	37	22.8	1	0.0
Head/Neck	5	3.1		
Bipolar Only			1	0.0
Cement Only			1	0.0
N Minor	53	32.7	6728	34.3
TOTAL	162	100.0	19643	100.0

Note: This table is restricted to revisions within 20.4 years for all groups to allow a time-matched comparison of revisions.

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 6**Revision Rates of Adapter (cless) Primary Total Conventional Hip Replacement by Fixation**

This analysis is provided as some prostheses have more than one fixation option. Additionally there are prostheses where an alternative to the recommended approach to fixation was used e.g. a cementless prosthesis that has been cemented or vice-versa.

Table 6: Revised Number of Adapter (cless) Primary Total Conventional Hip Replacement by Fixation

Fixation	N Revised	N Total
Cemented	0	1
Cementless	154	720
Hybrid (Femur Cemented)	8	16
Reverse Hybrid (Femur Cementless)	0	7
TOTAL	162	744

TABLE 7**Revision Rates of Adapter (cless) Primary Total Conventional Hip Replacement by Bearing Surface**

This analysis is provided as some prostheses are combined with a variety of bearing surfaces. All bearing surfaces used with this prosthesis are listed.

Table 7: Revised Number of Adapter (cless) Primary Total Conventional Hip Replacement by Bearing Surface

Bearing Surface	N Revised	N Total
Ceramic/Ceramic	10	50
Ceramic/Non XLPE	20	118
Ceramic/XLPE	2	39
Metal/Metal	92	363
Metal/Non XLPE	38	167
Metal/XLPE	0	7
TOTAL	162	744

TABLE 8**Revision Rates of Primary Total Conventional Hip Replacement by State**

This enables a state by state variation to be identified for the Adapter (cless) total conventional hip prosthesis and provides the comparative data for each of the states for all other total conventional hip prostheses.

The purpose of this analysis is to determine if the higher than anticipated rate of revision has widespread distribution between states. If there is widespread distribution then the reason for the higher than anticipated rate of revision is unlikely to be surgeon specific. If the prosthesis has been used in only a small number of states it is not possible to distinguish if the higher than anticipated rate of revision is related to the prosthesis, surgeon, technique or patient.

Table 8: Revised Number of Primary Total Conventional Hip Replacement by State

Component	State	N Revised	N Total
Adapter (cless)	NSW	84	464
	VIC	18	86
	QLD	7	18
	WA	16	41
	SA	1	6
	TAS	36	128
	ACT/NT	0	1
Other Total Conventional Hip	NSW	5507	167002
	VIC	4853	146558
	QLD	3968	104359
	WA	2510	64946
	SA	1887	52084
	TAS	439	18269
	ACT/NT	537	15394
TOTAL		19863	569356

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 9**Number of Revisions of Adapter (cless) Primary Total Conventional Hip Replacement by Year of Implant**

This analysis details the number of prostheses reported each year to the Registry for the Adapter (cless) total conventional hip prosthesis. It also provides the subsequent number of revisions of the primaries reported in that year.

Primary procedures performed in later years have had less follow up time therefore the number revised is expected to be less than the number revised in earlier years. For example, a primary procedure performed in 2025 has a maximum of one year to be revised, whereas a primary procedure performed in 2023 has a maximum of three years to be revised.

Table 9: Number of Revisions of Adapter (cless) Primary Total Conventional Hip Replacement by Year of Implant

Year of Implant	Number Revised	Total Number
2005	4	19
2006	36	140
2007	28	131
2008	27	122
2009	31	158
2010	30	115
2011	6	58
2013	0	1
TOTAL	162	744

TABLE 10**Revision Rates of Adapter (cless) Primary Total Conventional Hip Replacement by Catalogue Number Range**

Many prostheses have a number of catalogue ranges. The catalogue range is specific to particular design features; more than one catalogue range usually indicates a minor difference in design in a particular Adapter (cless) prosthesis.

This analysis has been undertaken to determine if the revision rate varies according to the catalogue number range.

Model	Catalogue Range	Catalogue Description	Cement	Material	Stem Type
Femoral Stem					
Adapter	11090302-11090373	COLLARLESS ADAPTER HIP STEM	NO	METAL	REQUIRES FEMNECK
Adapter	11099012-11099044	TINB-COATED COLLARLESS ADAPTER HIP STEM	NO	METAL	REQUIRES FEMNECK

Table 10: Revised Number of Adapter (cless) Primary Total Conventional Hip Replacement by Catalogue Number Range

Femoral Stem Range	N Revised	N Total
11090302-11090373	133	640
11099012-11099044	29	104
TOTAL	162	744

TABLE 11**Revision Rates of Adapter (cless) Primary Total Conventional Hip Replacement by Component**

A prosthesis may be combined with multiple components. This analysis has been undertaken to determine if the revision rate varies according to the component with which it is combined.

Table 11: Revised Number of Adapter (cless) Primary Total Conventional Hip Replacement by Acetabular Component

Acetabular Component	N Revised	N Total
2000	0	3
2000 Plus	8	79
Allofit	2	39
BSC-Cup	0	3
Bionik	138	527
Brunswick	0	2
DeltaLox	14	78
Eska Cup	0	4
PINNACLE	0	1
R3	0	1
Saturne	0	1
Trabecular Metal (Shell)	0	5
Trilogy	0	1
TOTAL	162	744