

Adapter (cemented) Total Conventional Hip Investigation

Note: This analysis compares the Adapter (ctd) femoral stem prosthesis with all other total conventional hip prostheses.

This prosthesis has been identified as having a significantly higher rate of revision. For a detailed explanation of the process used by the Registry that results in identification of prostheses that have a higher than anticipated rate of revision please refer to the Prostheses with Higher than Anticipated Rates of Revision chapter of the most recent AOANJRR Annual Report, <https://aoanjrr.sahmri.com/annual-reports-2026>.

Note: Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator. Procedures using prostheses with no recorded use in 2025 are excluded from the comparator.

TABLE 1

Revision Rate of Primary Total Conventional Hip Replacement

The revision rate of the Adapter (ctd) total conventional hip prosthesis is compared to all other total conventional hip prostheses.

Table 1: Revision Rates of Primary Total Conventional Hip Replacement

Component	N Revised	N Total	Obs. Years	Revisions/100 Obs. Yrs (95% CI)
Adapter (ctd)	34	148	1361	2.50 (1.73, 3.49)
Other Total Conventional Hip	19701	568612	3623071	0.54 (0.54, 0.55)
TOTAL	19735	568760	3624432	0.54 (0.54, 0.55)

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 2

Yearly Cumulative Percent Revision of Primary Total Conventional Hip Replacement

The yearly cumulative percent revision of the Adapter (ctd) total conventional hip prosthesis is compared to all other total conventional hip prostheses.

Table 2: Yearly Cumulative Percent Revision (95% CI) of Primary Total Conventional Hip Replacement

CPR	1 Yr	2 Yrs	3 Yrs	4 Yrs	5 Yrs	6 Yrs
Adapter (ctd)	4.1 (1.9, 8.9)	6.3 (3.3, 11.7)	9.2 (5.4, 15.3)	11.4 (7.2, 18.0)	16.3 (11.0, 23.7)	18.8 (13.1, 26.7)
Other Total Conventional Hip	1.7 (1.7, 1.8)	2.1 (2.1, 2.2)	2.5 (2.4, 2.5)	2.8 (2.7, 2.8)	3.0 (3.0, 3.1)	3.3 (3.2, 3.3)

CPR	7 Yrs	8 Yrs	9 Yrs	10 Yrs	11 Yrs	12 Yrs
Adapter (ctd)	21.6 (15.4, 29.8)	21.6 (15.4, 29.8)	22.6 (16.2, 31.0)	23.8 (17.1, 32.4)	25.1 (18.2, 34.0)	25.1 (18.2, 34.0)
Other Total Conventional Hip	3.5 (3.5, 3.6)	3.8 (3.7, 3.9)	4.1 (4.0, 4.2)	4.4 (4.3, 4.5)	4.7 (4.6, 4.8)	5.1 (5.0, 5.2)

CPR	13 Yrs	14 Yrs	15 Yrs	16 Yrs	17 Yrs	18 Yrs
Adapter (ctd)	26.7 (19.4, 36.0)	26.7 (19.4, 36.0)				
Other Total Conventional Hip	5.4 (5.3, 5.5)	5.8 (5.7, 5.9)	6.2 (6.1, 6.3)	6.6 (6.5, 6.7)	7.0 (6.8, 7.1)	7.3 (7.1, 7.5)

CPR	19 Yrs	20 Yrs	21 Yrs	22 Yrs	23 Yrs	24 Yrs	25 Yrs
Adapter (ctd)							
Other Total Conventional Hip	7.7 (7.5, 7.9)	8.0 (7.8, 8.2)	8.6 (8.3, 8.9)	9.0 (8.7, 9.4)	9.6 (9.1, 10.0)	9.6 (9.2, 10.1)	10.3 (9.4, 11.2)

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

FIGURE 1

Yearly Cumulative Percent Revision of Primary Total Conventional Hip Replacement

The yearly cumulative percent revision of the Adapter (ctd) total conventional hip prosthesis is compared to all other total conventional hip prostheses. In addition, hazard ratios are reported.

Hazard ratios are reported for specific time periods during which the hazard ratio is constant. This is done to enable more specific and valid comparisons of the risk of revision over time. The pattern of variation in risk has important implications with respect to the underlying reasons for any difference.

Figure 1: Cumulative Percent Revision of Primary Total Conventional Hip Replacement

***** Proportionality Not Achieved :-(*****

var1 VS var2 (Adapter (ctd) VS Other Total Conventional Hip) is non-proportional in time period 0_6Mth

The algorithm is not able to select a time point. Please specify time periods manually

The model will be fit over this entire period instead


var timepoint ProbChiSq

var1 0_6Mth 0.0279

var1 6Mth_max 0.4643

The PHREG Procedure

Dummy Text

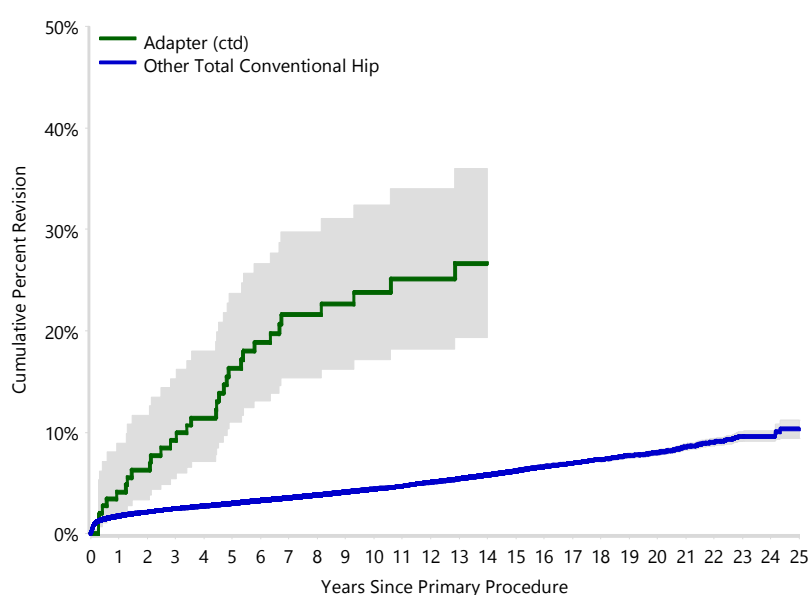

HR - adjusted for age and sex

Adapter (ctd) vs Other Total Conventional Hip

0 - 6Mth: HR=1.93 (95% CI 0.73, 5.15), p=0.187

6Mth+: HR=7.01 (95% CI 4.90, 10.03), p<0.001

Number at Risk	0 Yr	1 Yr	2 Yrs	3 Yrs	4 Yrs	5 Yrs	6 Yrs	7 Yrs
Adapter (ctd)	148	137	128	123	116	100	95	82
Other Total Conventional Hip	568612	501603	445433	391564	344389	299354	259516	220948

Number at Risk	8 Yrs	9 Yrs	10 Yrs	11 Yrs	12 Yrs	13 Yrs	14 Yrs	15 Yrs	16 Yrs
Adapter (ctd)	74	70	62	53	50	46	41	35	29
Other Total Conventional Hip	185798	153739	124848	99575	79802	63343	49730	38033	28355

Number at Risk	17 Yrs	18 Yrs	19 Yrs	20 Yrs	21 Yrs	22 Yrs	23 Yrs	24 Yrs	25 Yrs
Adapter (ctd)	24	15	4	1	0	0	0	0	0
Other Total Conventional Hip	20437	14638	10400	7396	4859	2863	1455	560	95

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 3**Primary Diagnosis for Revised Primary Total Conventional Hip Replacement**

This table identifies the diagnosis of the primary procedure which was subsequently revised. This information is provided as there is a variation on outcome depending on the primary diagnosis. It is therefore important when considering the reasons for a higher than anticipated rate of revision that there is identification of the primary diagnosis. This information should be compared to the primary diagnosis for the revisions of all other total conventional hip prostheses.

Table 3: Primary Diagnosis for Revised Primary Total Conventional Hip Replacement

Primary Diagnosis	Adapter (ctd)		Other Total Conventional Hip	
	Number	Percent	Number	Percent
Osteoarthritis	33	97.1	16296	82.7
Fractured Neck Of Femur			1521	7.7
Osteonecrosis			844	4.3
Developmental Dysplasia			321	1.6
Rheumatoid Arthritis			201	1.0
Failed Internal Fixation			165	0.8
Tumour			157	0.8
Other Inflammatory Arthritis	1	2.9	102	0.5
Fracture/Dislocation			60	0.3
Other			21	0.1
Arthrodesis Takedown			13	0.1
TOTAL	34	100.0	19701	100.0

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 4

Reasons for Revision

This is reported in two ways: a percentage of primary procedures revised and as a percentage of all revision procedures.

% Primaries Revised: This shows the proportional contribution of each revision diagnosis as a percentage of the total number of primary procedures. This percentage can be used to approximate the risk of being revised for that diagnosis. Differing percentages between groups, with the same distribution of follow up time, may identify problems of concern.

% Revisions: The number of revisions for each diagnosis is expressed as a percentage of the total number of revisions. This shows the distribution of reasons for revision within a group but cannot be used as a comparison between groups.

Table 4: Primary Total Conventional Hip Replacement - Reason for Revision (Follow-up Limited to 20.1 Years)

Revision Diagnosis	Number	Adapter (ctd)		Other Total Conventional Hip		
		% Primaries Revised	% Revisions	Number	% Primaries Revised	% Revisions
Infection	2	1.4	5.9	5008	0.9	25.5
Prosthesis Dislocation/Instability	8	5.4	23.5	4377	0.8	22.3
Fracture	1	0.7	2.9	4335	0.8	22.1
Loosening	15	10.1	44.1	3611	0.6	18.4
Pain	1	0.7	2.9	339	0.1	1.7
Leg Length Discrepancy				307	0.1	1.6
Malposition				282	0.0	1.4
Implant Breakage Stem				208	0.0	1.1
Lysis	1	0.7	2.9	196	0.0	1.0
Implant Breakage Acetabular Insert				126	0.0	0.6
Wear Acetabular Insert				106	0.0	0.5
Incorrect Sizing				96	0.0	0.5
Metal Related Pathology	6	4.1	17.6	90	0.0	0.5
Implant Breakage Acetabular				73	0.0	0.4
Wear Head				48	0.0	0.2
Tumour				46	0.0	0.2
Heterotopic Bone				26	0.0	0.1
Implant Breakage Head				26	0.0	0.1
Wear Acetabulum				6	0.0	0.0
Osteonecrosis				4	0.0	0.0
Synovitis				2	0.0	0.0
Progression Of Disease				1	0.0	0.0
Other				321	0.1	1.6
N Revision	34	23.0	100.0	19634	3.5	100.0
N Primary	148			568612		

Note: This table is restricted to revisions within 20.1 years for all groups to allow a time-matched comparison of revisions.

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

FIGURE 2**Cumulative Incidence Revision Diagnosis of Primary Total Conventional Hip Replacement**

This figure details the cumulative incidence of the most common reasons for revision. The five most common reasons for revision are included as long as each of these reasons account for more than 10 procedures or at least 5% of all revisions for the Adapter (ctd) total conventional hip prosthesis. A comparative graph is provided of the cumulative incidence for the same reasons for revisions for all other total conventional hip prostheses.

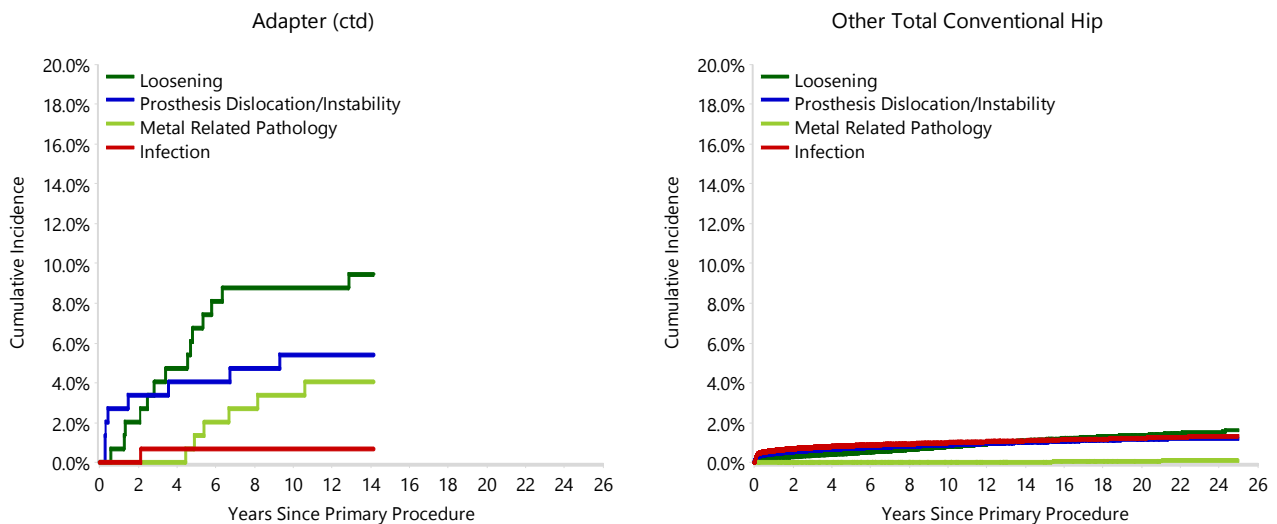
Figure 2: Cumulative Incidence Revision Diagnosis for Primary Total Conventional Hip Replacement

TABLE 5

Type of Revision Performed for Primary Total Conventional Hip Replacement

This analysis identifies the components used in the revision of the Adapter (ctd) total conventional hip prosthesis and compares it to the components used in the revision of all other total conventional hip prostheses.

The reason this analysis is undertaken is to identify whether there is one or more components which are being replaced that differ from the components replaced for revisions of all other total conventional hip prostheses i.e. is there a difference in the type of revision undertaken for the Adapter (ctd) total conventional hip prosthesis compared to all other total conventional hip prostheses.

Table 5: Primary Total Conventional Hip Replacement - Type of Revision (Follow-up Limited to 20.1 Years)

Type of Revision	Adapter (ctd)		Other Total Conventional Hip	
	Number	Percent	Number	Percent
Femoral Component	11	32.4	6630	33.8
Acetabular Component	9	26.5	3296	16.8
THR (Femoral/Acetabular)	5	14.7	2286	11.6
Cement Spacer	1	2.9	556	2.8
Removal of Prostheses			92	0.5
Reinsertion of Components			27	0.1
Total Femoral			12	0.1
Bipolar Head and Femoral			10	0.1
N Major	26	76.5	12909	65.7
Head/Insert	1	2.9	5357	27.3
Head Only			912	4.6
Minor Components			285	1.5
Insert Only			168	0.9
Head/Neck/Insert	6	17.6	1	0.0
Bipolar Only			1	0.0
Cement Only			1	0.0
Head/Neck	1	2.9		
N Minor	8	23.5	6725	34.3
TOTAL	34	100.0	19634	100.0

Note: This table is restricted to revisions within 20.1 years for all groups to allow a time-matched comparison of revisions.

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 6**Revision Rates of Adapter (ctd) Primary Total Conventional Hip Replacement by Fixation**

This analysis is provided as some prostheses have more than one fixation option. Additionally there are prostheses where an alternative to the recommended approach to fixation was used e.g. a cementless prosthesis that has been cemented or vice-versa.

Table 6: Revised Number of Adapter (ctd) Primary Total Conventional Hip Replacement by Fixation

Fixation	N Revised	N Total
Cemented	4	20
Cementless	0	1
Hybrid (Femur Cemented)	30	127
TOTAL	34	148

TABLE 7**Revision Rates of Adapter (ctd) Primary Total Conventional Hip Replacement by Bearing Surface**

This analysis is provided as some prostheses are combined with a variety of bearing surfaces. All bearing surfaces used with this prosthesis are listed.

Table 7: Revised Number of Adapter (ctd) Primary Total Conventional Hip Replacement by Bearing Surface

Bearing Surface	N Revised	N Total
Ceramic/Ceramic	7	24
Ceramic/Non XLPE	3	23
Metal/Metal	18	59
Metal/Non XLPE	6	39
Metal/XLPE	0	3
TOTAL	34	148

TABLE 8

Revision Rates of Primary Total Conventional Hip Replacement by State

This enables a state by state variation to be identified for the Adapter (ctd) total conventional hip prosthesis and provides the comparative data for each of the states for all other total conventional hip prostheses.

The purpose of this analysis is to determine if the higher than anticipated rate of revision has widespread distribution between states. If there is widespread distribution then the reason for the higher than anticipated rate of revision is unlikely to be surgeon specific. If the prosthesis has been used in only a small number of states it is not possible to distinguish if the higher than anticipated rate of revision is related to the prosthesis, surgeon, technique or patient.

Table 8: Revised Number of Primary Total Conventional Hip Replacement by State

Component	State	N Revised	N Total
Adapter (ctd)	NSW	23	98
	VIC	0	8
	QLD	0	2
	TAS	11	40
Other Total Conventional Hip	NSW	5507	167002
	VIC	4853	146558
	QLD	3968	104359
	WA	2510	64946
	SA	1887	52084
	TAS	439	18269
	ACT/NT	537	15394
TOTAL		19735	568760

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 9**Number of Revisions of Adapter (ctd) Primary Total Conventional Hip Replacement by Year of Implant**

This analysis details the number of prostheses reported each year to the Registry for the Adapter (ctd) total conventional hip prosthesis. It also provides the subsequent number of revisions of the primaries reported in that year.

Primary procedures performed in later years have had less follow up time therefore the number revised is expected to be less than the number revised in earlier years. For example, a primary procedure performed in 2025 has a maximum of one year to be revised, whereas a primary procedure performed in 2023 has a maximum of three years to be revised.

Table 9: Number of Revisions of Adapter (ctd) Primary Total Conventional Hip Replacement by Year of Implant

Year of Implant	Number Revised	Total Number
2005	3	7
2006	9	41
2007	14	52
2008	6	33
2009	2	8
2010	0	7
TOTAL	34	148

TABLE 10**Revision Rates of Adapter (ctd) Primary Total Conventional Hip Replacement by Catalogue Number Range**

Many prostheses have a number of catalogue ranges. The catalogue range is specific to particular design features; more than one catalogue range usually indicates a minor difference in design in a particular Adapter (ctd) prosthesis.

This analysis has been undertaken to determine if the revision rate varies according to the catalogue number range.

Model	Catalogue Range	Catalogue Description	Cement	Material	Stem Type	Fixation
Femoral Stem						
Adapter	11030002-11030043	CEMENTED COLLARLESS ADAPTER HIP STEM	YES	METAL	REQUIRES FEMNECK	MATT

Table 10: Revised Number of Adapter (ctd) Primary Total Conventional Hip Replacement by Catalogue Number Range

Femoral Stem Range	N Revised	N Total
11030002-11030043	34	148
TOTAL	34	148

TABLE 11**Revision Rates of Adapter (ctd) Primary Total Conventional Hip Replacement by Component**

A prosthesis may be combined with multiple components. This analysis has been undertaken to determine if the revision rate varies according to the component with which it is combined.

Table 11: Revised Number of Adapter (ctd) Primary Total Conventional Hip Replacement by Acetabular Component

Acetabular Component	N Revised	N Total
2000	0	1
2000 Plus	11	53
Allofit	0	2
Bionik	20	72
Eska Cup	3	19
Reflection (Shell)	0	1
TOTAL	34	148