

Optetrak Logic CR/Optetrak Logic Total Knee Investigation

Note: This analysis compares the Optetrak Logic CR/Optetrak Logic femoral/tibial combination with all other total knee prostheses.

This combination has been identified as having a significantly higher rate of revision. For a detailed explanation of the process used by the Registry that results in identification of prostheses that have a higher than anticipated rate of revision please refer to the Prostheses with Higher than Anticipated Rates of Revision chapter of the most recent AOANJRR Annual Report, <https://aoanjrr.sahmri.com/annual-reports-2026>.

Note: Procedures using prostheses with no recorded use in 2025 are excluded from the comparator.

TABLE 1

Revision Rate of Primary Total Knee Replacement

The revision rate of the Optetrak Logic CR/Optetrak Logic total knee combination is compared to all other total knee prostheses.

Table 1: Revision Rates of Primary Total Knee Replacement

Component	N Revised	N Total	Obs. Years	Revisions/100 Obs. Yrs (95% CI)
Optetrak Logic CR/Optetrak Logic	256	1887	12471	2.05 (1.81, 2.32)
Other Total Knee	27284	791430	5145922	0.53 (0.52, 0.54)
TOTAL	27540	793317	5158393	0.53 (0.53, 0.54)

Note: Prostheses no longer used in 2025 are excluded from the comparator.

TABLE 2

Yearly Cumulative Percent Revision of Primary Total Knee Replacement

The yearly cumulative percent revision of the Optetrak Logic CR/Optetrak Logic total knee combination is compared to all other total knee prostheses.

Table 2: Yearly Cumulative Percent Revision (95% CI) of Primary Total Knee Replacement

CPR	1 Yr	2 Yrs	3 Yrs	4 Yrs	5 Yrs	6 Yrs
Optetrak Logic CR/Optetrak Logic	1.1 (0.7, 1.6)	2.3 (1.7, 3.1)	4.3 (3.5, 5.4)	6.3 (5.2, 7.5)	8.3 (7.1, 9.7)	10.3 (9.0, 11.9)
Other Total Knee	1.0 (1.0, 1.0)	1.8 (1.8, 1.8)	2.4 (2.3, 2.4)	2.8 (2.7, 2.8)	3.1 (3.1, 3.1)	3.4 (3.4, 3.5)

CPR	7 Yrs	8 Yrs	9 Yrs	10 Yrs	11 Yrs	12 Yrs
Optetrak Logic CR/Optetrak Logic	13.4 (11.8, 15.2)	15.8 (13.9, 17.9)	18.6 (16.4, 21.2)	21.5 (18.6, 24.7)		
Other Total Knee	3.7 (3.7, 3.8)	4.0 (4.0, 4.1)	4.3 (4.3, 4.4)	4.6 (4.6, 4.7)	4.9 (4.9, 5.0)	5.3 (5.2, 5.3)

CPR	13 Yrs	14 Yrs	15 Yrs	16 Yrs	17 Yrs	18 Yrs
Optetrak Logic CR/Optetrak Logic						
Other Total Knee	5.6 (5.5, 5.7)	5.9 (5.8, 6.0)	6.3 (6.2, 6.4)	6.7 (6.6, 6.9)	7.2 (7.0, 7.3)	7.5 (7.4, 7.7)

CPR	19 Yrs	20 Yrs	21 Yrs	22 Yrs	23 Yrs	24 Yrs	25 Yrs
Optetrak Logic CR/Optetrak Logic							
Other Total Knee	7.9 (7.7, 8.0)	8.1 (7.9, 8.3)	8.5 (8.2, 8.7)	8.7 (8.4, 9.0)	8.8 (8.6, 9.2)	9.1 (8.7, 9.5)	9.5 (8.8, 10.2)

Note: Prostheses no longer used in 2025 are excluded from the comparator.


FIGURE 1

Yearly Cumulative Percent Revision of Primary Total Knee Replacement

The yearly cumulative percent revision of the Optetrak Logic CR/Optetrak Logic total knee combination is compared to all other total knee prostheses. In addition, hazard ratios are reported.

Hazard ratios are reported for specific time periods during which the hazard ratio is constant. This is done to enable more specific and valid comparisons of the risk of revision over time. The pattern of variation in risk has important implications with respect to the underlying reasons for any difference.

Figure 1: Cumulative Percent Revision of Primary Total Knee Replacement

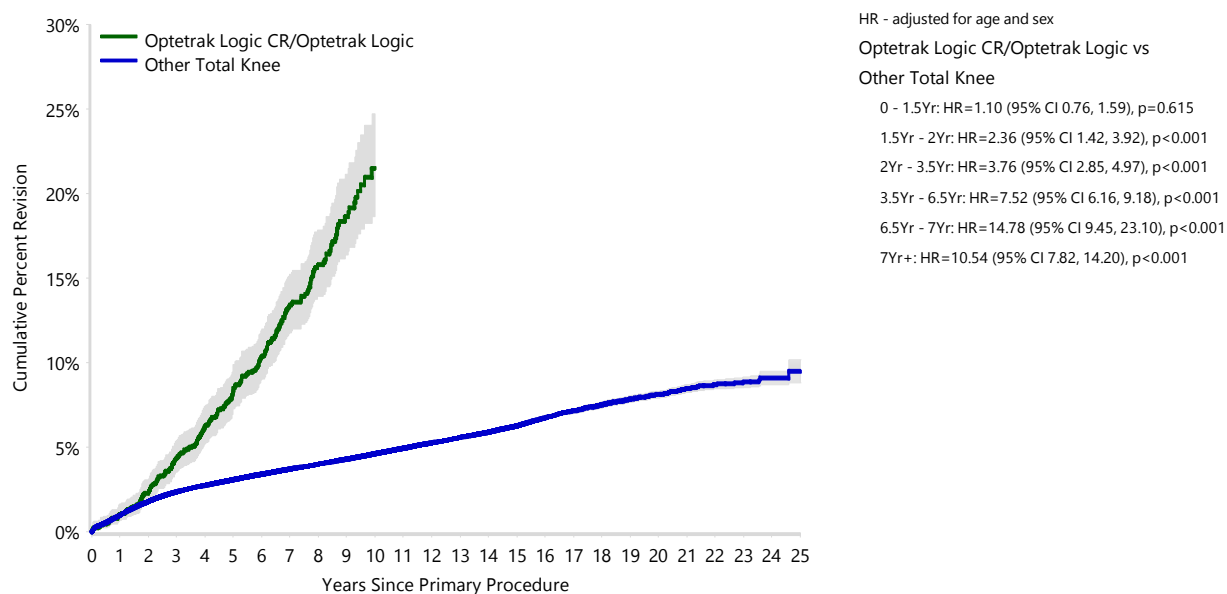

Number at Risk	0 Yr	1 Yr	2 Yrs	3 Yrs	4 Yrs	5 Yrs	6 Yrs	7 Yrs
Optetrak Logic CR/Optetrak Logic	1887	1861	1826	1750	1623	1369	1162	891
Other Total Knee	791430	708699	628274	550835	487544	424868	369845	315999

Number at Risk	8 Yrs	9 Yrs	10 Yrs	11 Yrs	12 Yrs	13 Yrs	14 Yrs	15 Yrs	16 Yrs
Optetrak Logic CR/Optetrak Logic	580	321	134	15	0	0	0	0	0
Other Total Knee	266816	222320	182654	146980	117201	92187	70380	52526	38351

Number at Risk	17 Yrs	18 Yrs	19 Yrs	20 Yrs	21 Yrs	22 Yrs	23 Yrs	24 Yrs	25 Yrs
Optetrak Logic CR/Optetrak Logic	0	0	0	0	0	0	0	0	0
Other Total Knee	27628	19434	13616	9427	6164	3657	2096	972	244

Note: Prostheses no longer used in 2025 are excluded from the comparator.

TABLE 3**Primary Diagnosis for Revised Primary Total Knee Replacement**

This table identifies the diagnosis of the primary procedure which was subsequently revised. This information is provided as there is a variation on outcome depending on the primary diagnosis. It is therefore important when considering the reasons for a higher than anticipated rate of revision that there is identification of the primary diagnosis. This information should be compared to the primary diagnosis for the revisions of all other total knee prostheses.

Table 3: Primary Diagnosis for Revised Primary Total Knee Replacement

Primary Diagnosis	Optetrak Logic CR/Optetrak Logic		Other Total Knee	
	Number	Percent	Number	Percent
Osteoarthritis	248	96.9	26392	96.7
Rheumatoid Arthritis	3	1.2	330	1.2
Tumour			199	0.7
Other Inflammatory Arthritis	4	1.6	177	0.6
Osteonecrosis	1	0.4	104	0.4
Fracture			55	0.2
Other			26	0.1
Chondrocalcinosis			1	0.0
TOTAL	256	100.0	27284	100.0

Note: Prostheses no longer used in 2025 are excluded from the comparator.

TABLE 4

Reasons for Revision

This is reported in two ways: a percentage of primary procedures revised and as a percentage of all revision procedures.

% Primaries Revised: This shows the proportional contribution of each revision diagnosis as a percentage of the total number of primary procedures. This percentage can be used to approximate the risk of being revised for that diagnosis. Differing percentages between groups, with the same distribution of follow up time, may identify problems of concern.

% Revisions: The number of revisions for each diagnosis is expressed as a percentage of the total number of revisions. This shows the distribution of reasons for revision within a group but cannot be used as a comparison between groups.

Table 4: Primary Total Knee Replacement - Reason for Revision (Follow-up Limited to 11.2 Years)

Revision Diagnosis	Optetrak Logic CR/Optetrak Logic			Other Total Knee		
	Number	% Primaries Revised	% Revisions	Number	% Primaries Revised	% Revisions
Infection	26	1.4	10.2	7612	1.0	29.9
Loosening	70	3.7	27.3	5318	0.7	20.9
Instability	17	0.9	6.6	2738	0.3	10.8
Patella Erosion	4	0.2	1.6	1785	0.2	7.0
Pain	5	0.3	2.0	1766	0.2	6.9
Patellofemoral Pain				1473	0.2	5.8
Arthrofibrosis	3	0.2	1.2	1168	0.1	4.6
Fracture	1	0.1	0.4	988	0.1	3.9
Malalignment	3	0.2	1.2	503	0.1	2.0
Lysis	54	2.9	21.1	233	0.0	0.9
Wear Tibial Insert	40	2.1	15.6	231	0.0	0.9
Incorrect Sizing	1	0.1	0.4	217	0.0	0.9
Patella Maltracking	1	0.1	0.4	191	0.0	0.8
Implant Breakage Tibial Insert	16	0.8	6.3	167	0.0	0.7
Bearing Dislocation				142	0.0	0.6
Implant Breakage Patella				140	0.0	0.6
Metal Related Pathology	2	0.1	0.8	106	0.0	0.4
Prosthesis Dislocation				75	0.0	0.3
Synovitis	5	0.3	2.0	62	0.0	0.2
Osteonecrosis				50	0.0	0.2
Wear Patella	2	0.1	0.8	41	0.0	0.2
Implant Breakage Femoral	1	0.1	0.4	36	0.0	0.1
Implant Breakage Tibial				36	0.0	0.1
Tumour				33	0.0	0.1
Heterotopic Bone				14	0.0	0.1
Progression Of Disease				7	0.0	0.0
Wear Tibial				5	0.0	0.0
Patella Dislocation				3	0.0	0.0
Incorrect Side				1	0.0	0.0
Wear Femoral				1	0.0	0.0
Other	5	0.3	2.0	308	0.0	1.2
N Revision	256	13.6	100.0	25450	3.2	100.0
N Primary	1887			791430		

Note: This table is restricted to revisions within 11.2 years for all groups to allow a time-matched comparison of revisions.

Note: Prostheses no longer used in 2025 are excluded from the comparator.

FIGURE 2**Cumulative Incidence Revision Diagnosis of Primary Total Knee Replacement**

This figure details the cumulative incidence of the most common reasons for revision. The five most common reasons for revision are included as long as each of these reasons account for more than 10 procedures or at least 5% of all revisions for the Optetrak Logic CR/Optetrak Logic total knee combination. A comparative graph is provided of the cumulative incidence for the same reasons for revisions for all other total knee prostheses.

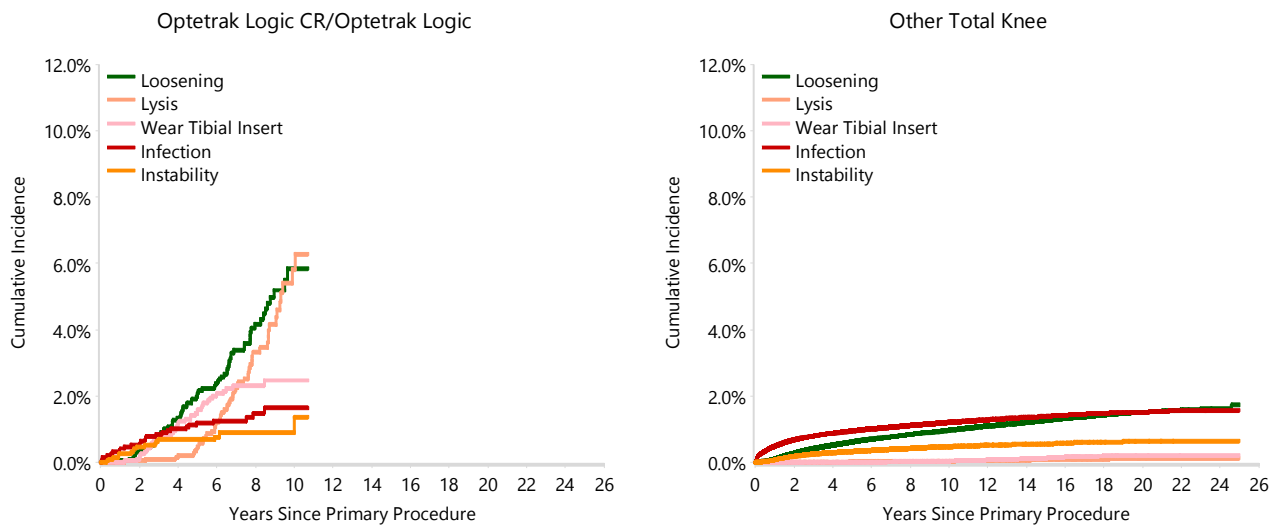
Figure 2: Cumulative Incidence Revision Diagnosis for Primary Total Knee Replacement

TABLE 5

Type of Revision Performed for Primary Total Knee Replacement

This analysis identifies the components used in the revision of the Optetrak Logic CR/Optetrak Logic total knee combination and compares it to the components used in the revision of all other total knee prostheses.

The reason this analysis is undertaken is to identify whether there is one or more components which are being replaced that differ from the components replaced for revisions of all other total knee prostheses i.e. is there a difference in the type of revision undertaken for the Optetrak Logic CR/Optetrak Logic total knee combination compared to all other total knee prostheses.

Table 5: Primary Total Knee Replacement - Type of Revision (Follow-up Limited to 11.2 Years)

Type of Revision	Optetrak Logic CR/Optetrak Logic		Other Total Knee	
	Number	Percent	Number	Percent
TKR (Tibial/Femoral)	111	43.4	6022	23.7
Tibial Component	1	0.4	1859	7.3
Femoral Component	35	13.7	1275	5.0
Cement Spacer	2	0.8	1103	4.3
Removal of Prostheses			130	0.5
Total Femoral			20	0.1
Reinsertion of Components			8	0.0
N Major	149	58.2	10417	40.9
Insert Only	73	28.5	8217	32.3
Patella Only	7	2.7	4293	16.9
Insert/Patella	27	10.5	2450	9.6
Minor Components			64	0.3
Cement Only			9	0.0
N Minor	107	41.8	15033	59.1
TOTAL	256	100.0	25450	100.0

Note: This table is restricted to revisions within 11.2 years for all groups to allow a time-matched comparison of revisions.

Note: Prostheses no longer used in 2025 are excluded from the comparator.

TABLE 6**Revision Rates of Optetrak Logic CR/Optetrak Logic Primary Total Knee Replacement by Fixation**

This analysis is provided as some prostheses have more than one fixation option. Additionally there are prostheses where an alternative to the recommended approach to fixation was used e.g. a cementless prosthesis that has been cemented or vice-versa.

Table 6: Revised Number of Optetrak Logic CR/Optetrak Logic Primary Total Knee Replacement by Fixation

Fixation	N Revised	N Total
Cemented	74	758
Hybrid (Tibial Cemented)	182	1129
TOTAL	256	1887

TABLE 7**Revision Rates of Optetrak Logic CR/Optetrak Logic Primary Total Knee Replacement by Bearing Surface**

This analysis is provided as some prostheses are combined with a variety of bearing surfaces. All bearing surfaces used with this combination are listed.

Table 7: Revised Number of Optetrak Logic CR/Optetrak Logic Primary Total Knee Replacement by Bearing Surface

Bearing Surface	N Revised	N Total
Non XLPE	256	1887
TOTAL	256	1887

TABLE 8**Revision Rates of Optetrak Logic CR/Optetrak Logic Primary Total Knee Replacement by Bearing Mobility**

This analysis is provided as some prostheses are combined with a variety of bearing mobilities. All bearing mobilities used with this combination are listed.

Table 8: Revised Number of Optetrak Logic CR/Optetrak Logic Primary Total Knee Replacement by Bearing Mobility

Bearing Mobility	N Revised	N Total
Fixed	256	1887
TOTAL	256	1887

TABLE 9**Revision Rates of Optetrak Logic CR/Optetrak Logic Primary Total Knee Replacement by Stability**

This analysis is provided as some prostheses are combined with a variety of stabilities. All stabilities used with this combination are listed.

Table 9: Revised Number of Optetrak Logic CR/Optetrak Logic Primary Total Knee Replacement by Stability

Stability	N Revised	N Total
Minimally Stabilised	256	1887
TOTAL	256	1887

TABLE 10

Revision Rates of Primary Total Knee Replacement by State

This enables a state by state variation to be identified for the Optetrak Logic CR/Optetrak Logic total knee combination and provides the comparative data for each of the states for all other total knee prostheses.

The purpose of this analysis is to determine if the higher than anticipated rate of revision has widespread distribution between states. If there is widespread distribution then the reason for the higher than anticipated rate of revision is unlikely to be surgeon specific. If the prosthesis has been used in only a small number of states it is not possible to distinguish if the higher than anticipated rate of revision is related to the prosthesis, surgeon, technique or patient.

Table 10: Revised Number of Primary Total Knee Replacement by State

Component	State	N Revised	N Total
Optetrak Logic CR/Optetrak Logic	NSW	191	1366
	VIC	11	83
	QLD	36	316
	WA	16	85
	TAS	2	37
Other Total Knee	NSW	7550	259497
	VIC	6020	163565
	QLD	5789	167244
	WA	3600	92841
	SA	3220	73114
	TAS	475	15075
	ACT/NT	630	20094
TOTAL		27540	793317

Note: Prostheses no longer used in 2025 are excluded from the comparator.

TABLE 11**Number of Revisions of Optetrak Logic CR/Optetrak Logic Primary Total Knee Replacement by Year of Implant**

This analysis details the number of prostheses reported each year to the Registry for the Optetrak Logic CR/Optetrak Logic total knee combination. It also provides the subsequent number of revisions of the primaries reported in that year.

Primary procedures performed in later years have had less follow up time therefore the number revised is expected to be less than the number revised in earlier years. For example, a primary procedure performed in 2025 has a maximum of one year to be revised, whereas a primary procedure performed in 2023 has a maximum of three years to be revised.

Table 11: Number of Revisions of Optetrak Logic CR/Optetrak Logic Primary Total Knee Replacement by Year of Implant

Year of Implant	Number Revised	Total Number
2014	0	21
2015	19	171
2016	36	234
2017	45	321
2018	39	339
2019	55	286
2020	27	195
2021	28	226
2022	7	83
2023	0	11
TOTAL	256	1887

TABLE 12

Revision Rates of Optetrak Logic CR/Optetrak Logic Primary Total Knee Replacement by Catalogue Number Range

Many prostheses have a number of catalogue ranges. The catalogue range is specific to particular design features; more than one catalogue range usually indicates a minor difference in design in a particular Optetrak Logic CR/Optetrak Logic prosthesis.

This analysis has been undertaken to determine if the revision rate varies according to the catalogue number range.

Model	Catalogue Range	Catalogue Description	Cement
Femoral			
Optetrak Logic CR	02010030200-02010030360	OPTETRAK LOGIC FEMORAL COMPONENT CR CEMENTED	YES
Optetrak Logic CR	02010040200-02010040360	OPTETRAK LOGIC FEMORAL COMPONENT CR POROUS	NO
Tibial			
Optetrak Logic	02012450000-02012456060	OPTETRAK LOGIC FIT TIBIAL TRAY CEMENTED	YES

Table 12: Revised Number of Optetrak Logic CR/Optetrak Logic Primary Total Knee Replacement by Catalogue Number Range

Femoral Range	Tibial Range	N Revised	N Total
02010030200-02010030360	02012450000-02012456060	74	752
02010040200-02010040360	02012450000-02012456060	182	1135
TOTAL		256	1887