

BI-MENTUM (cemented cup) Total Conventional Hip Investigation

Note: This analysis compares the BI-MENTUM (cemented cup) acetabular prosthesis with all other total conventional hip prostheses.

This prosthesis has been identified as having a significantly higher rate of revision. For a detailed explanation of the process used by the Registry that results in identification of prostheses that have a higher than anticipated rate of revision please refer to the Prostheses with Higher than Anticipated Rates of Revision chapter of the most recent AOANJRR Annual Report, <https://aoanjrr.sahmri.com/annual-reports-2026>.

Note: Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator. Procedures using prostheses with no recorded use in 2025 are excluded from the comparator.

TABLE 1

Revision Rate of Primary Total Conventional Hip Replacement

The revision rate of the BI-MENTUM (cemented cup) total conventional hip prosthesis is compared to all other total conventional hip prostheses.

Table 1: Revision Rates of Primary Total Conventional Hip Replacement

Component	N Revised	N Total	Obs. Years	Revisions/100 Obs. Yrs (95% CI)
BI-MENTUM (cemented cup)	11	75	128	8.58 (4.28, 15.35)
Other Total Conventional Hip	19690	568537	3622943	0.54 (0.54, 0.55)
TOTAL	19701	568612	3623071	0.54 (0.54, 0.55)

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 2

Yearly Cumulative Percent Revision of Primary Total Conventional Hip Replacement

The yearly cumulative percent revision of the BI-MENTUM (cemented cup) total conventional hip prosthesis is compared to all other total conventional hip prostheses.

Table 2: Yearly Cumulative Percent Revision (95% CI) of Primary Total Conventional Hip Replacement

CPR	1 Yr	2 Yrs	3 Yrs	4 Yrs	5 Yrs	6 Yrs
BI-MENTUM (cemented cup)	14.6 (8.1, 25.5)	17.1 (9.7, 29.1)	17.1 (9.7, 29.1)			
Other Total Conventional Hip	1.7 (1.7, 1.8)	2.1 (2.1, 2.2)	2.5 (2.4, 2.5)	2.7 (2.7, 2.8)	3.0 (3.0, 3.1)	3.3 (3.2, 3.3)

CPR	7 Yrs	8 Yrs	9 Yrs	10 Yrs	11 Yrs	12 Yrs
BI-MENTUM (cemented cup)						
Other Total Conventional Hip	3.5 (3.5, 3.6)	3.8 (3.7, 3.9)	4.1 (4.0, 4.2)	4.4 (4.3, 4.5)	4.7 (4.6, 4.8)	5.1 (5.0, 5.2)

CPR	13 Yrs	14 Yrs	15 Yrs	16 Yrs	17 Yrs	18 Yrs
BI-MENTUM (cemented cup)						
Other Total Conventional Hip	5.4 (5.3, 5.5)	5.8 (5.7, 5.9)	6.2 (6.1, 6.3)	6.6 (6.5, 6.7)	7.0 (6.8, 7.1)	7.3 (7.1, 7.5)

CPR	19 Yrs	20 Yrs	21 Yrs	22 Yrs	23 Yrs	24 Yrs	25 Yrs
BI-MENTUM (cemented cup)							
Other Total Conventional Hip	7.7 (7.5, 7.9)	8.0 (7.8, 8.2)	8.6 (8.3, 8.9)	9.0 (8.7, 9.4)	9.6 (9.1, 10.0)	9.6 (9.2, 10.1)	10.3 (9.4, 11.2)

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.


FIGURE 1

Yearly Cumulative Percent Revision of Primary Total Conventional Hip Replacement

The yearly cumulative percent revision of the BI-MENTUM (cemented cup) total conventional hip prosthesis is compared to all other total conventional hip prostheses. In addition, hazard ratios are reported.

Hazard ratios are reported for specific time periods during which the hazard ratio is constant. This is done to enable more specific and valid comparisons of the risk of revision over time. The pattern of variation in risk has important implications with respect to the underlying reasons for any difference.

Figure 1: Cumulative Percent Revision of Primary Total Conventional Hip Replacement

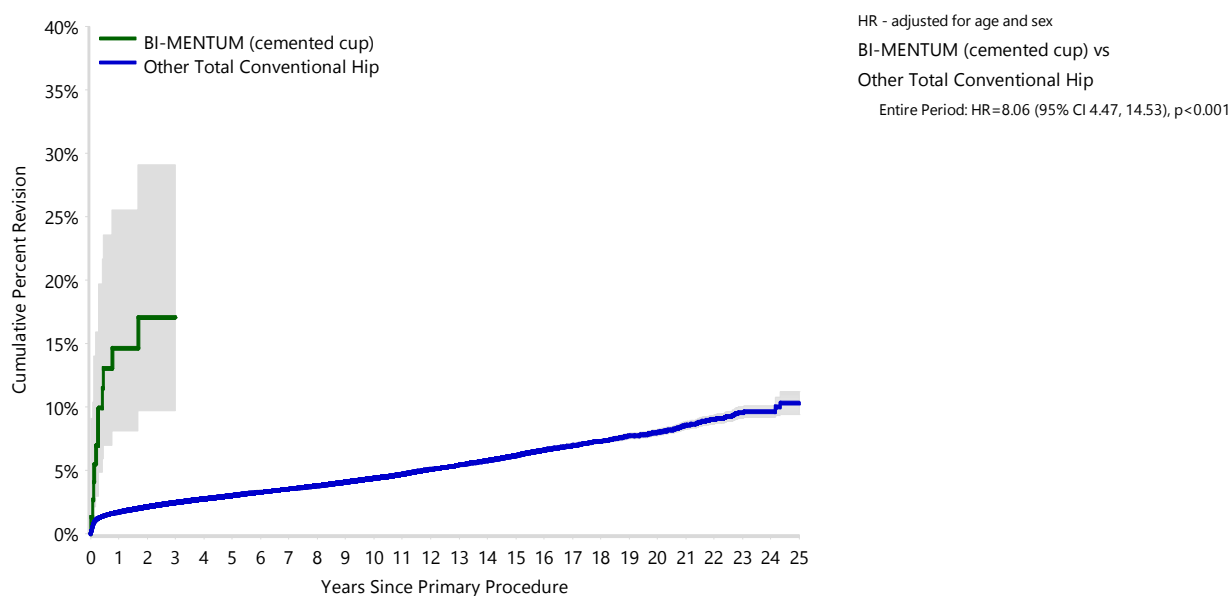

Number at Risk	0 Yr	1 Yr	2 Yrs	3 Yrs	4 Yrs	5 Yrs	6 Yrs	7 Yrs
BI-MENTUM (cemented cup)	75	49	29	12	6	0	0	0
Other Total Conventional Hip	568537	501554	445404	391552	344383	299354	259516	220948

Number at Risk	8 Yrs	9 Yrs	10 Yrs	11 Yrs	12 Yrs	13 Yrs	14 Yrs	15 Yrs	16 Yrs
BI-MENTUM (cemented cup)	0	0	0	0	0	0	0	0	0
Other Total Conventional Hip	185798	153739	124848	99575	79802	63343	49730	38033	28355

Number at Risk	17 Yrs	18 Yrs	19 Yrs	20 Yrs	21 Yrs	22 Yrs	23 Yrs	24 Yrs	25 Yrs
BI-MENTUM (cemented cup)	0	0	0	0	0	0	0	0	0
Other Total Conventional Hip	20437	14638	10400	7396	4859	2863	1455	560	95

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 3**Primary Diagnosis for Revised Primary Total Conventional Hip Replacement**

This table identifies the diagnosis of the primary procedure which was subsequently revised. This information is provided as there is a variation on outcome depending on the primary diagnosis. It is therefore important when considering the reasons for a higher than anticipated rate of revision that there is identification of the primary diagnosis. This information should be compared to the primary diagnosis for the revisions of all other total conventional hip prostheses.

Table 3: Primary Diagnosis for Revised Primary Total Conventional Hip Replacement

Primary Diagnosis	BI-MENTUM (cemented cup)		Other Total Conventional Hip	
	Number	Percent	Number	Percent
Osteoarthritis	5	45.5	16291	82.7
Fractured Neck Of Femur	1	9.1	1520	7.7
Osteonecrosis	3	27.3	841	4.3
Developmental Dysplasia			321	1.6
Rheumatoid Arthritis			201	1.0
Failed Internal Fixation			165	0.8
Tumour			157	0.8
Other Inflammatory Arthritis	1	9.1	101	0.5
Fracture/Dislocation	1	9.1	59	0.3
Other			21	0.1
Arthrodesis Takedown			13	0.1
TOTAL	11	100.0	19690	100.0

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 4

Reasons for Revision

This is reported in two ways: a percentage of primary procedures revised and as a percentage of all revision procedures.

% Primaries Revised: This shows the proportional contribution of each revision diagnosis as a percentage of the total number of primary procedures. This percentage can be used to approximate the risk of being revised for that diagnosis. Differing percentages between groups, with the same distribution of follow up time, may identify problems of concern.

% Revisions: The number of revisions for each diagnosis is expressed as a percentage of the total number of revisions. This shows the distribution of reasons for revision within a group but cannot be used as a comparison between groups.

Table 4: Primary Total Conventional Hip Replacement - Reason for Revision (Follow-up Limited to 4.5 Years)

Revision Diagnosis	BI-MENTUM (cemented cup)			Other Total Conventional Hip		
	Number	% Primaries Revised	% Revisions	Number	% Primaries Revised	% Revisions
Infection	7	9.3	63.6	4353	0.8	30.1
Prosthesis Dislocation/Instability				3652	0.6	25.2
Fracture	2	2.7	18.2	2811	0.5	19.4
Loosening	2	2.7	18.2	2198	0.4	15.2
Leg Length Discrepancy				284	0.0	2.0
Pain				253	0.0	1.7
Malposition				224	0.0	1.5
Incorrect Sizing				89	0.0	0.6
Implant Breakage Acetabular Insert				86	0.0	0.6
Implant Breakage Stem				68	0.0	0.5
Implant Breakage Acetabular				45	0.0	0.3
Lysis				45	0.0	0.3
Tumour				34	0.0	0.2
Metal Related Pathology				26	0.0	0.2
Heterotopic Bone				22	0.0	0.2
Wear Head				17	0.0	0.1
Wear Acetabular Insert				12	0.0	0.1
Implant Breakage Head				9	0.0	0.1
Osteonecrosis				4	0.0	0.0
Wear Acetabulum				2	0.0	0.0
Progression Of Disease				1	0.0	0.0
Synovitis				1	0.0	0.0
Other				239	0.0	1.7
N Revision	11	14.7	100.0	14475	2.5	100.0
N Primary	75			568537		

Note: This table is restricted to revisions within 4.5 years for all groups to allow a time-matched comparison of revisions.

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

FIGURE 2**Cumulative Incidence Revision Diagnosis of Primary Total Conventional Hip Replacement**

This figure details the cumulative incidence of the most common reasons for revision. The five most common reasons for revision are included as long as each of these reasons account for more than 10 procedures or at least 5% of all revisions for the BI-MENTUM (cemented cup) total conventional hip prosthesis. A comparative graph is provided of the cumulative incidence for the same reasons for revisions for all other total conventional hip prostheses.

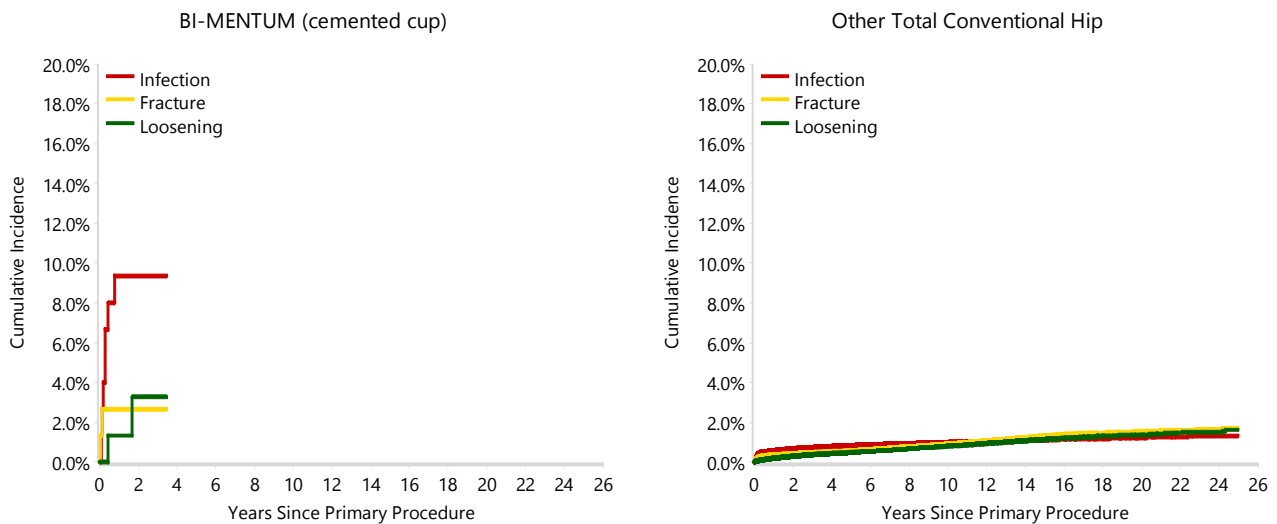
Figure 2: Cumulative Incidence Revision Diagnosis for Primary Total Conventional Hip Replacement

TABLE 5

Type of Revision Performed for Primary Total Conventional Hip Replacement

This analysis identifies the components used in the revision of the BI-MENTUM (cemented cup) total conventional hip prosthesis and compares it to the components used in the revision of all other total conventional hip prostheses.

The reason this analysis is undertaken is to identify whether there is one or more components which are being replaced that differ from the components replaced for revisions of all other total conventional hip prostheses i.e. is there a difference in the type of revision undertaken for the BI-MENTUM (cemented cup) total conventional hip prosthesis compared to all other total conventional hip prostheses.

Table 5: Primary Total Conventional Hip Replacement - Type of Revision (Follow-up Limited to 4.5 Years)

Type of Revision	BI-MENTUM (cemented cup)		Other Total Conventional Hip	
	Number	Percent	Number	Percent
Femoral Component			4279	29.6
Acetabular Component	1	9.1	2427	16.8
THR (Femoral/Acetabular)	8	72.7	1534	10.6
Cement Spacer			431	3.0
Removal of Prostheses			78	0.5
Reinsertion of Components			26	0.2
Bipolar Head and Femoral			8	0.1
Total Femoral			8	0.1
N Major	9	81.8	8791	60.7
Head/Insert	1	9.1	4521	31.2
Head Only			808	5.6
Minor Components	1	9.1	199	1.4
Insert Only			153	1.1
Bipolar Only			1	0.0
Cement Only			1	0.0
Head/Neck/Insert			1	0.0
N Minor	2	18.2	5684	39.3
TOTAL	11	100.0	14475	100.0

Note: This table is restricted to revisions within 4.5 years for all groups to allow a time-matched comparison of revisions.

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 6**Revision Rates of BI-MENTUM (cemented cup) Primary Total Conventional Hip Replacement by Fixation**

This analysis is provided as some prostheses have more than one fixation option. Additionally there are prostheses where an alternative to the recommended approach to fixation was used e.g. a cementless prosthesis that has been cemented or vice-versa.

Table 6: Revised Number of BI-MENTUM (cemented cup) Primary Total Conventional Hip Replacement by Fixation

Fixation	N Revised	N Total
Cemented	10	57
Reverse Hybrid (Femur Cementless)	1	18
TOTAL	11	75

TABLE 7**Revision Rates of BI-MENTUM (cemented cup) Primary Total Conventional Hip Replacement by Bearing Surface**

This analysis is provided as some prostheses are combined with a variety of bearing surfaces. All bearing surfaces used with this prosthesis are listed.

Table 7: Revised Number of BI-MENTUM (cemented cup) Primary Total Conventional Hip Replacement by Bearing Surface

Bearing Surface	N Revised	N Total
Ceramic/Non XLPE	2	25
Metal/Non XLPE	9	50
TOTAL	11	75

TABLE 8**Revision Rates of BI-MENTUM (cemented cup) Primary Total Conventional Hip Replacement by Approach**

This analysis is provided as some prostheses are used with a variety of surgical approaches. All surgical approaches used with this prosthesis are listed.

Table 8: Revised Number of BI-MENTUM (cemented cup) Primary Total Conventional Hip Replacement by Approach

Approach	N Revised	N Total
Anterior	0	3
Lateral	0	7
Posterior	11	64
TOTAL	11	74

Note: Excludes 1 procedures with no approach recorded

TABLE 9

Revision Rates of Primary Total Conventional Hip Replacement by State

This enables a state by state variation to be identified for the BI-MENTUM (cemented cup) total conventional hip prosthesis and provides the comparative data for each of the states for all other total conventional hip prostheses.

The purpose of this analysis is to determine if the higher than anticipated rate of revision has widespread distribution between states. If there is widespread distribution then the reason for the higher than anticipated rate of revision is unlikely to be surgeon specific. If the prosthesis has been used in only a small number of states it is not possible to distinguish if the higher than anticipated rate of revision is related to the prosthesis, surgeon, technique or patient.

Table 9: Revised Number of Primary Total Conventional Hip Replacement by State

Component	State	N Revised	N Total
BI-MENTUM (cemented cup)	NSW	3	27
	VIC	0	3
	QLD	0	6
	WA	7	34
	SA	1	5
Other Total Conventional Hip	NSW	5504	166975
	VIC	4853	146555
	QLD	3968	104353
	WA	2503	64912
	SA	1886	52079
	TAS	439	18269
	ACT/NT	537	15394
TOTAL		19701	568612

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 10**Number of Revisions of BI-MENTUM (cemented cup) Primary Total Conventional Hip Replacement by Year of Implant**

This analysis details the number of prostheses reported each year to the Registry for the BI-MENTUM (cemented cup) total conventional hip prosthesis. It also provides the subsequent number of revisions of the primaries reported in that year.

Primary procedures performed in later years have had less follow up time therefore the number revised is expected to be less than the number revised in earlier years. For example, a primary procedure performed in 2025 has a maximum of one year to be revised, whereas a primary procedure performed in 2023 has a maximum of three years to be revised.

Table 10: Number of Revisions of BI-MENTUM (cemented cup) Primary Total Conventional Hip Replacement by Year of Implant

Year of Implant	Number Revised	Total Number
2021	1	8
2022	2	14
2023	3	23
2024	2	21
2025	3	9
TOTAL	11	75

TABLE 11

Revision Rates of BI-MENTUM (cemented cup) Primary Total Conventional Hip Replacement by Catalogue Number Range

Many prostheses have a number of catalogue ranges. The catalogue range is specific to particular design features; more than one catalogue range usually indicates a minor difference in design in a particular BI-MENTUM (cemented cup) prosthesis.

This analysis has been undertaken to determine if the revision rate varies according to the catalogue number range.

Model	Catalogue Range	Catalogue Description	Cement	Material
Acetabular				
BI-MENTUM	DS4901043-DS4901063	CEMENTED ACETABULAR CUP	YES	METAL

Table 11: Revised Number of BI-MENTUM (cemented cup) Primary Total Conventional Hip Replacement by Catalogue Number Range

Acetabular Range	N Revised	N Total
DS4901043-DS4901063	11	75
TOTAL	11	75

TABLE 12

Revision Rates of BI-MENTUM (cemented cup) Primary Total Conventional Hip Replacement by Component

A prosthesis may be combined with multiple components. This analysis has been undertaken to determine if the revision rate varies according to the component with which it is combined.

Table 12: Revised Number of BI-MENTUM (cemented cup) Primary Total Conventional Hip Replacement by Femoral Stem Component

Femoral Stem Component	N Revised	N Total
ACTIS	0	7
Arcos	0	1
C-Stem AMT	5	38
CORAIL	2	15
CPCS	0	1
Exeter V40	2	7
LPS	0	1
Reclaim	0	1
Short Exeter V40	2	3
VerSys	0	1
TOTAL	11	75