

Arcos (cone body) Total Conventional Hip Investigation

Note: This analysis compares the Arcos (cone body) femoral stem prosthesis with all other total conventional hip prostheses.

This prosthesis has been identified as having a significantly higher rate of revision. For a detailed explanation of the process used by the Registry that results in identification of prostheses that have a higher than anticipated rate of revision please refer to the Prostheses with Higher than Anticipated Rates of Revision chapter of the most recent AOANJRR Annual Report, <https://aoanjrr.sahmri.com/annual-reports-2026>.

Note: Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator. Procedures using prostheses with no recorded use in 2025 are excluded from the comparator.

TABLE 1

Revision Rate of Primary Total Conventional Hip Replacement

The revision rate of the Arcos (cone body) total conventional hip prosthesis is compared to all other total conventional hip prostheses.

Table 1: Revision Rates of Primary Total Conventional Hip Replacement

Component	N Revised	N Total	Obs. Years	Revisions/100 Obs. Yrs (95% CI)
Arcos (cone body)	34	413	1412	2.41 (1.67, 3.36)
Other Total Conventional Hip	19672	568281	3622188	0.54 (0.54, 0.55)
TOTAL	19706	568694	3623600	0.54 (0.54, 0.55)

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 2

Yearly Cumulative Percent Revision of Primary Total Conventional Hip Replacement

The yearly cumulative percent revision of the Arcos (cone body) total conventional hip prosthesis is compared to all other total conventional hip prostheses.

Table 2: Yearly Cumulative Percent Revision (95% CI) of Primary Total Conventional Hip Replacement

CPR	1 Yr	2 Yrs	3 Yrs	4 Yrs	5 Yrs	6 Yrs
Arcos (cone body)	6.4 (4.4, 9.4)	8.1 (5.7, 11.4)	9.1 (6.4, 12.7)	9.1 (6.4, 12.7)	9.1 (6.4, 12.7)	10.1 (7.0, 14.4)
Other Total Conventional Hip	1.7 (1.7, 1.8)	2.1 (2.1, 2.2)	2.5 (2.4, 2.5)	2.7 (2.7, 2.8)	3.0 (3.0, 3.1)	3.3 (3.2, 3.3)

CPR	7 Yrs	8 Yrs	9 Yrs	10 Yrs	11 Yrs	12 Yrs
Arcos (cone body)	10.1 (7.0, 14.4)	10.1 (7.0, 14.4)				
Other Total Conventional Hip	3.5 (3.5, 3.6)	3.8 (3.7, 3.9)	4.1 (4.0, 4.2)	4.4 (4.3, 4.5)	4.7 (4.6, 4.8)	5.1 (5.0, 5.2)

CPR	13 Yrs	14 Yrs	15 Yrs	16 Yrs	17 Yrs	18 Yrs
Arcos (cone body)						
Other Total Conventional Hip	5.4 (5.3, 5.5)	5.8 (5.7, 5.9)	6.2 (6.1, 6.3)	6.6 (6.5, 6.7)	7.0 (6.8, 7.1)	7.3 (7.1, 7.5)

CPR	19 Yrs	20 Yrs	21 Yrs	22 Yrs	23 Yrs	24 Yrs	25 Yrs
Arcos (cone body)							
Other Total Conventional Hip	7.7 (7.5, 7.9)	8.0 (7.8, 8.2)	8.6 (8.3, 8.8)	9.0 (8.7, 9.4)	9.6 (9.1, 10.0)	9.6 (9.1, 10.1)	10.3 (9.4, 11.2)

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.


FIGURE 1

Yearly Cumulative Percent Revision of Primary Total Conventional Hip Replacement

The yearly cumulative percent revision of the Arcos (cone body) total conventional hip prosthesis is compared to all other total conventional hip prostheses. In addition, hazard ratios are reported.

Hazard ratios are reported for specific time periods during which the hazard ratio is constant. This is done to enable more specific and valid comparisons of the risk of revision over time. The pattern of variation in risk has important implications with respect to the underlying reasons for any difference.

Figure 1: Cumulative Percent Revision of Primary Total Conventional Hip Replacement

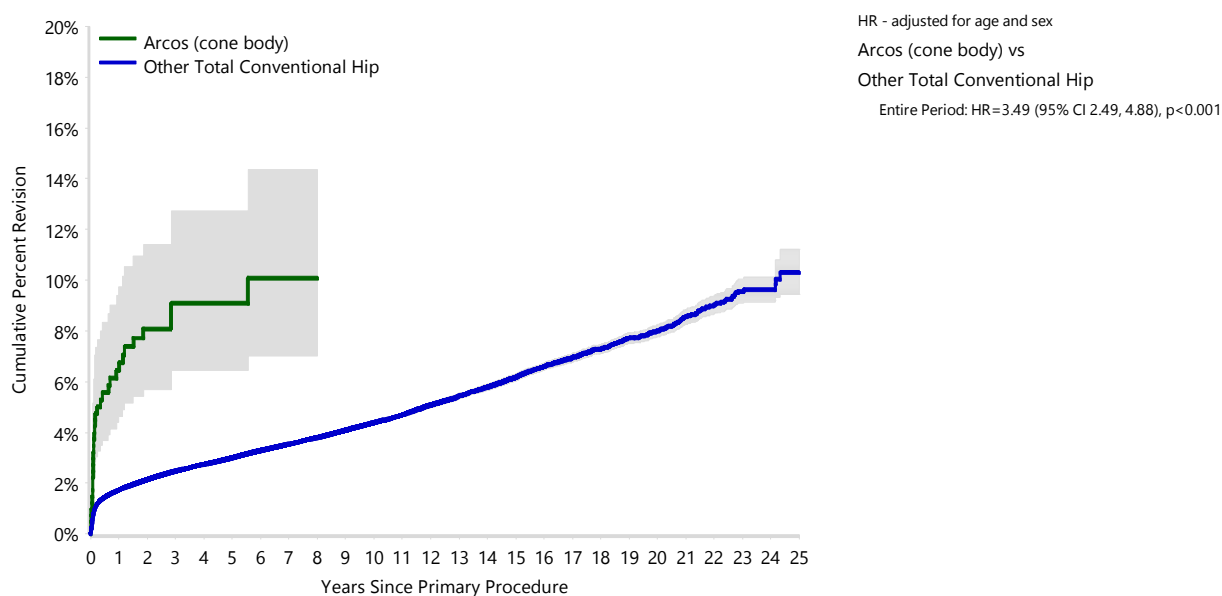

Number at Risk	0 Yr	1 Yr	2 Yrs	3 Yrs	4 Yrs	5 Yrs	6 Yrs	7 Yrs
Arcos (cone body)	413	308	233	176	132	105	79	65
Other Total Conventional Hip	568281	501367	445265	391447	344311	299298	259482	220925

Number at Risk	8 Yrs	9 Yrs	10 Yrs	11 Yrs	12 Yrs	13 Yrs	14 Yrs	15 Yrs	16 Yrs
Arcos (cone body)	44	27	19	17	12	6	3	0	0
Other Total Conventional Hip	185787	153734	124845	99574	79801	63342	49729	38033	28355

Number at Risk	17 Yrs	18 Yrs	19 Yrs	20 Yrs	21 Yrs	22 Yrs	23 Yrs	24 Yrs	25 Yrs
Arcos (cone body)	0	0	0	0	0	0	0	0	0
Other Total Conventional Hip	20437	14638	10400	7396	4859	2863	1455	560	95

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 3**Primary Diagnosis for Revised Primary Total Conventional Hip Replacement**

This table identifies the diagnosis of the primary procedure which was subsequently revised. This information is provided as there is a variation on outcome depending on the primary diagnosis. It is therefore important when considering the reasons for a higher than anticipated rate of revision that there is identification of the primary diagnosis. This information should be compared to the primary diagnosis for the revisions of all other total conventional hip prostheses.

Table 3: Primary Diagnosis for Revised Primary Total Conventional Hip Replacement

Primary Diagnosis	Arcos (cone body)		Other Total Conventional Hip	
	Number	Percent	Number	Percent
Osteoarthritis	5	14.7	16291	82.8
Fractured Neck Of Femur	9	26.5	1513	7.7
Osteonecrosis	1	2.9	843	4.3
Developmental Dysplasia	2	5.9	319	1.6
Rheumatoid Arthritis			201	1.0
Tumour	1	2.9	156	0.8
Failed Internal Fixation	14	41.2	155	0.8
Other Inflammatory Arthritis			102	0.5
Fracture/Dislocation	1	2.9	59	0.3
Other			21	0.1
Arthrodesis Takedown	1	2.9	12	0.1
TOTAL	34	100.0	19672	100.0

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 4

Reasons for Revision

This is reported in two ways: a percentage of primary procedures revised and as a percentage of all revision procedures.

% Primaries Revised: This shows the proportional contribution of each revision diagnosis as a percentage of the total number of primary procedures. This percentage can be used to approximate the risk of being revised for that diagnosis. Differing percentages between groups, with the same distribution of follow up time, may identify problems of concern.

% Revisions: The number of revisions for each diagnosis is expressed as a percentage of the total number of revisions. This shows the distribution of reasons for revision within a group but cannot be used as a comparison between groups.

Table 4: Primary Total Conventional Hip Replacement - Reason for Revision (Follow-up Limited to 14.3 Years)

Revision Diagnosis	Arcos (cone body)			Other Total Conventional Hip		
	Number	% Primaries Revised	% Revisions	Number	% Primaries Revised	% Revisions
Infection	13	3.1	38.2	4936	0.9	25.9
Prosthesis Dislocation/Instability	11	2.7	32.4	4300	0.8	22.5
Fracture	2	0.5	5.9	4183	0.7	21.9
Loosening	3	0.7	8.8	3481	0.6	18.2
Pain	1	0.2	2.9	331	0.1	1.7
Leg Length Discrepancy				307	0.1	1.6
Malposition				280	0.0	1.5
Implant Breakage Stem	1	0.2	2.9	194	0.0	1.0
Lysis				168	0.0	0.9
Implant Breakage Acetabular Insert				123	0.0	0.6
Incorrect Sizing				96	0.0	0.5
Wear Acetabular Insert				81	0.0	0.4
Metal Related Pathology	1	0.2	2.9	73	0.0	0.4
Implant Breakage Acetabular				70	0.0	0.4
Tumour				46	0.0	0.2
Wear Head				44	0.0	0.2
Heterotopic Bone	1	0.2	2.9	26	0.0	0.1
Implant Breakage Head				23	0.0	0.1
Wear Acetabulum				6	0.0	0.0
Osteonecrosis				4	0.0	0.0
Synovitis				2	0.0	0.0
Progression Of Disease				1	0.0	0.0
Other	1	0.2	2.9	318	0.1	1.7
N Revision	34	8.2	100.0	19093	3.4	100.0
N Primary	413			568281		

Note: This table is restricted to revisions within 14.3 years for all groups to allow a time-matched comparison of revisions.

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

FIGURE 2**Cumulative Incidence Revision Diagnosis of Primary Total Conventional Hip Replacement**

This figure details the cumulative incidence of the most common reasons for revision. The five most common reasons for revision are included as long as each of these reasons account for more than 10 procedures or at least 5% of all revisions for the Arcos (cone body) total conventional hip prosthesis. A comparative graph is provided of the cumulative incidence for the same reasons for revisions for all other total conventional hip prostheses.

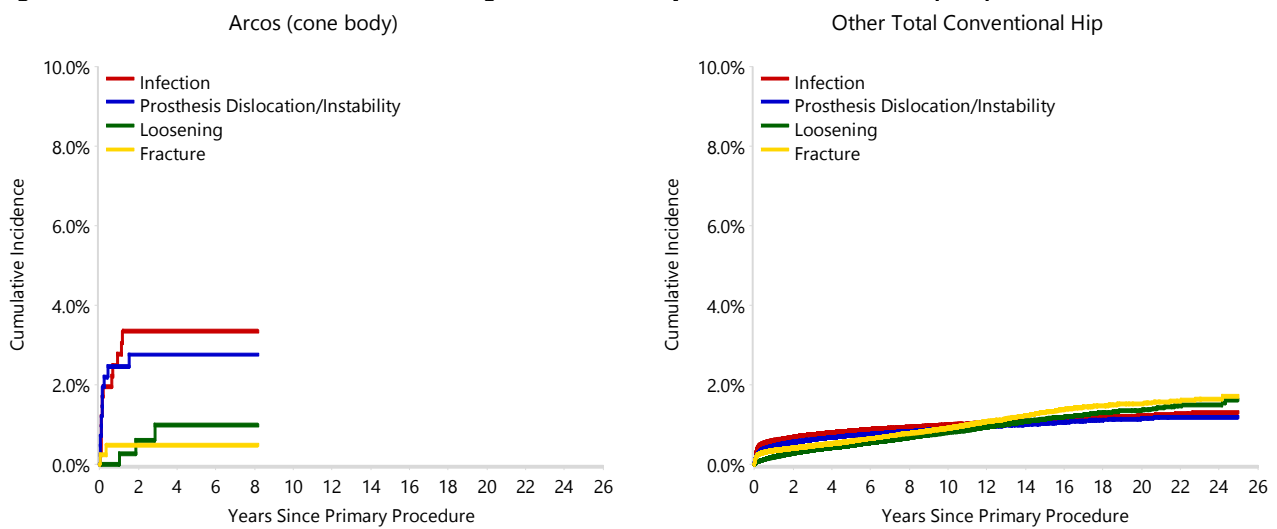
Figure 2: Cumulative Incidence Revision Diagnosis for Primary Total Conventional Hip Replacement

TABLE 5

Type of Revision Performed for Primary Total Conventional Hip Replacement

This analysis identifies the components used in the revision of the Arcos (cone body) total conventional hip prosthesis and compares it to the components used in the revision of all other total conventional hip prostheses.

The reason this analysis is undertaken is to identify whether there is one or more components which are being replaced that differ from the components replaced for revisions of all other total conventional hip prostheses i.e. is there a difference in the type of revision undertaken for the Arcos (cone body) total conventional hip prosthesis compared to all other total conventional hip prostheses.

Table 5: Primary Total Conventional Hip Replacement - Type of Revision (Follow-up Limited to 14.3 Years)

Type of Revision	Arcos (cone body)		Other Total Conventional Hip	
	Number	Percent	Number	Percent
Femoral Component	10	29.4	6422	33.6
Acetabular Component	4	11.8	3185	16.7
THR (Femoral/Acetabular)	6	17.6	2189	11.5
Cement Spacer	1	2.9	548	2.9
Removal of Prostheses			90	0.5
Reinsertion of Components			27	0.1
Total Femoral			12	0.1
Bipolar Head and Femoral			10	0.1
N Major	21	61.8	12483	65.4
Head/Insert	13	38.2	5262	27.6
Head Only			903	4.7
Minor Components			277	1.5
Insert Only			165	0.9
Bipolar Only			1	0.0
Cement Only			1	0.0
Head/Neck/Insert			1	0.0
N Minor	13	38.2	6610	34.6
TOTAL	34	100.0	19093	100.0

Note: This table is restricted to revisions within 14.3 years for all groups to allow a time-matched comparison of revisions.

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 6**Revision Rates of Arcos (cone body) Primary Total Conventional Hip Replacement by Fixation**

This analysis is provided as some prostheses have more than one fixation option. Additionally there are prostheses where an alternative to the recommended approach to fixation was used e.g. a cementless prosthesis that has been cemented or vice-versa.

Table 6: Revised Number of Arcos (cone body) Primary Total Conventional Hip Replacement by Fixation

Fixation	N Revised	N Total
Cemented	0	2
Cementless	34	380
Hybrid (Femur Cemented)	0	2
Reverse Hybrid (Femur Cementless)	0	29
TOTAL	34	413

TABLE 7**Revision Rates of Arcos (cone body) Primary Total Conventional Hip Replacement by Bearing Surface**

This analysis is provided as some prostheses are combined with a variety of bearing surfaces. All bearing surfaces used with this prosthesis are listed.

Table 7: Revised Number of Arcos (cone body) Primary Total Conventional Hip Replacement by Bearing Surface

Bearing Surface	N Revised	N Total
Ceramic/Ceramic	1	8
Ceramic/Non XLPE	0	1
Ceramic/XLPE	1	15
Ceramic/XLPE + Antioxidant	12	88
Metal/Non XLPE	1	41
Metal/XLPE	4	24
Metal/XLPE + Antioxidant	15	234
Unknown	0	2
TOTAL	34	413

TABLE 8**Revision Rates of Arcos (cone body) Primary Total Conventional Hip Replacement by Approach**

This analysis is provided as some prostheses are used with a variety of surgical approaches. All surgical approaches used with this prosthesis are listed.

Table 8: Revised Number of Arcos (cone body) Primary Total Conventional Hip Replacement by Approach

Approach	N Revised	N Total
Anterior	0	12
Lateral	3	49
Posterior	26	299
TOTAL	29	360

Note: Excludes 53 procedures with no approach recorded

TABLE 9

Revision Rates of Primary Total Conventional Hip Replacement by State

This enables a state by state variation to be identified for the Arcos (cone body) total conventional hip prosthesis and provides the comparative data for each of the states for all other total conventional hip prostheses.

The purpose of this analysis is to determine if the higher than anticipated rate of revision has widespread distribution between states. If there is widespread distribution then the reason for the higher than anticipated rate of revision is unlikely to be surgeon specific. If the prosthesis has been used in only a small number of states it is not possible to distinguish if the higher than anticipated rate of revision is related to the prosthesis, surgeon, technique or patient.

Table 9: Revised Number of Primary Total Conventional Hip Replacement by State

Component	State	N Revised	N Total
Arcos (cone body)	NSW	6	86
	VIC	12	96
	QLD	7	144
	WA	3	54
	SA	4	26
	TAS	2	5
	ACT/NT	0	2
Other Total Conventional Hip	NSW	5502	166934
	VIC	4843	146473
	QLD	3963	104264
	WA	2507	64892
	SA	1883	52062
	TAS	437	18264
	ACT/NT	537	15392
TOTAL		19706	568694

Note: Prostheses no longer used in 2025 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

TABLE 10**Number of Revisions of Arcos (cone body) Primary Total Conventional Hip Replacement by Year of Implant**

This analysis details the number of prostheses reported each year to the Registry for the Arcos (cone body) total conventional hip prosthesis. It also provides the subsequent number of revisions of the primaries reported in that year.

Primary procedures performed in later years have had less follow up time therefore the number revised is expected to be less than the number revised in earlier years. For example, a primary procedure performed in 2025 has a maximum of one year to be revised, whereas a primary procedure performed in 2023 has a maximum of three years to be revised.

Table 10: Number of Revisions of Arcos (cone body) Primary Total Conventional Hip Replacement by Year of Implant

Year of Implant	Number Revised	Total Number
2011	1	9
2012	1	14
2013	2	14
2014	1	12
2015	0	9
2016	3	15
2017	1	26
2018	2	34
2019	1	17
2020	2	32
2021	2	25
2022	4	45
2023	2	46
2024	5	59
2025	7	56
TOTAL	34	413

TABLE 11

Revision Rates of Arcos (cone body) Primary Total Conventional Hip Replacement by Catalogue Number Range

Many prostheses have a number of catalogue ranges. The catalogue range is specific to particular design features; more than one catalogue range usually indicates a minor difference in design in a particular Arcos (cone body) prosthesis.

This analysis has been undertaken to determine if the revision rate varies according to the catalogue number range.

Model	Catalogue Range	Catalogue Description	Cement	Material
Femoral Stem				
Arcos	11301300-11301310	CONE PROXIMAL BODY 50MM	NO	METAL
Arcos	11301301-11301347	CONE PROXIMAL BODY W/STANDARD OFFSET >50MM	NO	METAL
Arcos	11301311-11301357	CONE PROXIMAL BODY W/HIGH OFFSET >50MM	NO	METAL
Arcos	22301300-22301345	CONE PROXIMAL BODY TYPE 1 W/STANDARD OFFSET	NO	METAL
Arcos	22301310-22301355	CONE PROXIMAL BODY TYPE 1 W/HIGH OFFSET	NO	METAL

Table 11: Revised Number of Arcos (cone body) Primary Total Conventional Hip Replacement by Catalogue Number Range

Femoral Stem Range	N Revised	N Total
11301300-11301310	6	106
11301301-11301347	11	193
11301311-11301357	10	61
22301300-22301345	5	45
22301310-22301355	2	8
TOTAL	34	413

TABLE 12

Revision Rates of Arcos (cone body) Primary Total Conventional Hip Replacement by Component

A prosthesis may be combined with multiple components. This analysis has been undertaken to determine if the revision rate varies according to the component with which it is combined.

Table 12: Revised Number of Arcos (cone body) Primary Total Conventional Hip Replacement by Acetabular Component

Acetabular Component	N Revised	N Total
Acetabular Shell (Global)	1	2
Allofit	0	1
Avantage	1	13
BI-MENTUM (ctd)	0	1
Continuum	1	14
Custom Made (Ossis)	0	1
Delta Revision TT	0	1
Delta-TT	0	1
Exceed	0	9
Freedom	0	9
G7	24	274
Jump System	0	1
Logical G	0	1
Low Profile Cup	0	1
Mallory-Head	0	6
Mpact	0	2
Muller	0	1
Novae	0	1
PINNACLE	0	2
R3	1	1
Recap	0	1
Regenerex	2	46
Trabecular Metal (Shell)	2	12
Tri-Flange	0	1
Trident (Shell)	0	1
Trident II/Tritanium (Shell) (dual mobility)	0	1
Trident/Tritanium (Shell)	1	1
Trilogy IT	1	1
Trinity	0	5
Versafitcup CC	0	1
ZCA	0	1
TOTAL	34	413